(12) INTERNATIONAL APPLICATION PUBLISHED UNDER THE PATENT COOPERATION TREATY (PCT)

(19) World Intellectual Property Organization International Bureau





(43) International Publication Date 30 May 2002 (30.05.2002)

(51) International Patent Classification7:

PCT

C12N 15/12,

(10) International Publication Number WO~02/42461~A2

C07K 14/705

(21) International Application Number: PCT/US01/44386

(22) International Filing Date:

26 November 2001 (26.11.2001)

(25) Filing Language: English

(26) Publication Language: English


(30) Priority Data: 60/253,404 27 November 2000 (27.11.2000) US 60/255,366 12 December 2000 (12.12.2000) US 20 February 2001 (20.02.2001) 60/270,266 US 60/270,286 20 February 2001 (20.02.2001) US 60/282,032 6 April 2001 (06.04.2001) US 60/282,356 6 April 2001 (06.04.2001) US 60/282,358 6 April 2001 (06.04.2001) US 60/282,365 6 April 2001 (06.04.2001) US 60/290,917 14 May 2001 (14.05.2001) US 31 July 2001 (31.07.2001) 60/309,208 US

- (71) Applicant (for all designated States except US): ARENA PHARMACEUTICALS, INC. [US/US]; 6166 Nancy Drive, San Diego, CA 92121 (US).
- (72) Inventors; and
- (75) Inventors/Applicants (for US only): CHEN, Ruoping [CN/US]; 5296 Timber Branch Way, San Diego, CA 92150 (US). CHU, Zhi, Liang [CN/US]; 11489 Creekstone Lane, San Diego, CA 92128 (US). DANG, Huong, T. [VN/US]; 5352 Oak Park Drive, San Diego, CA 92105 (US). LOWITZ, Kevin, P. [US/US]; 8031 Caminito de Pizza, Apt. C, San Diego, CA 92108 (US). PRIDE, Cameron [US/US]; 11026 West Ocean Air Drive #1128, San Diego, CA 92130 (US).
- (74) Agents: STRAHER, Michael, P. et al.; Woodcock Washburn LLP, 46th floor, One Liberty Place, Philadelphia, PA 19103 (US).
- (81) Designated States (national): AE, AG, AL, AM, AT, AU, AZ, BA, BB, BG, BR, BY, BZ, CA, CH, CN, CO, CR, CU, CZ, DE, DK, DM, DZ, EC, EE, ES, FI, GB, GD, GE, GH, GM, HR, HU, ID, IL, IN, IS, JP, KE, KG, KP, KR, KZ, LC, LK, LR, LS, LT, LU, LV, MA, MD, MG, MK, MN, MW,

[Continued on next page]

(54) Title: ENDOGENOUS AND NON-ENDOGENOUS VERSIONS OF HUMAN G PROTEIN-COUPLED RECEPTORS

IP3 Assay in 293 Cells

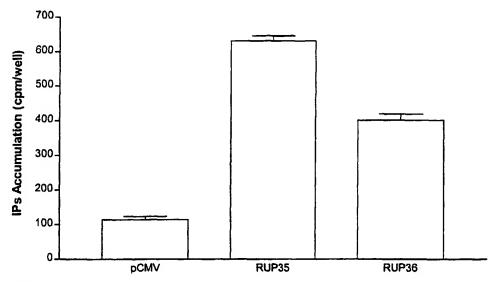

(57) Abstract: The invention disclosed in this patent document relates to transmembrane receptors, more particularly to a human G protein-coupled receptor for which the endogenous ligand is unknown, and to mutated (non-endogenous) versions of the human GPCRs for evidence of constitutive activity.



02/42461 A2

MX, MZ, NO, NZ, PL, PT, RO, RU, SD, SE, SG, SI, SK, SL, TJ, TM, TR, TT, TZ, UA, UG, US, UZ, VN, YU, ZA, ZW.

(84) Designated States (regional): ARIPO patent (GH, GM, KE, LS, MW, MZ, SD, SL, SZ, TZ, UG, ZM, ZW), Eurasian patent (AM, AZ, BY, KG, KZ, MD, RU, TJ, TM), European patent (AT, BE, CH, CY, DE, DK, ES, FI, FR, GB, GR, IE, IT, LU, MC, NL, PT, SE, TR), OAPI patent (BF, BJ, CF, CG, CI, CM, GA, GN, GQ, GW, ML, MR, NE, SN, TD, TG).

Declaration under Rule 4.17:

- of inventorship (Rule 4.17(iv)) for US only

Published:

 without international search report and to be republished upon receipt of that report

For two-letter codes and other abbreviations, refer to the "Guidance Notes on Codes and Abbreviations" appearing at the beginning of each regular issue of the PCT Gazette.

ENDOGENOUS AND NON-ENDOGENOUS VERSIONS OF HUMAN G PROTEIN-COUPLED RECEPTORS

CROSS-REFERENCE TO RELATED APPLICATIONS

This application is a continuation-in-part of U.S. Serial Number 09/170,496, filed with the United States Patent and Trademark Office on October 13, 1998 and its corresponding PCT application number PCT/US99/23938, published as WO 00/22129 on April 20, 2000. This document claims the benefit of priority from the following provisional applications, all filed via U.S. Express Mail with the United States Patent and Trademark Office on the indicated dates: U.S. Provisional Number 60/253,404, filed November 27, 2000; U.S. Provisional Number 60/255,366, filed December 12, 2000; U.S. Provisional Number 60/270,286 filed February 20, 2001; U.S. Provisional Number 60/282,356, filed April 6, 2001, which claims priority from U.S. Provisional Number 60/270,266, filed February 20, 2001; U.S. Provisional Number 60/282,032, filed April 6, 2001; U.S. Provisional Number 60/282,355, filed April 6, 2001; U.S. Provisional Number 60/290,917, filed May 14, 2001; U.S. Provisional Number 60/309,208, filed July 31, 2001; the disclosures of which are incorporated in their entirety by reference.

FIELD OF THE INVENTION

5

The present invention relates to transmembrane receptors, in some embodiments to G protein-coupled receptors and, in some preferred embodiments, to endogenous GPCRs that are altered to establish or enhance constitutive activity of the receptor. In some embodiments, the constitutively activated GPCRs will be used for the direct identification of candidate compounds as receptor agonists or inverse agonists having applicability as therapeutic agents.

BACKGROUND OF THE INVENTION

Although a number of receptor classes exist in humans, by far the most abundant and therapeutically relevant is represented by the G protein-coupled receptor (GPCR) class. It is estimated that there are some 30,000-40,000 genes within the human genome, and of these, approximately 2% are estimated to code for GPCRs. Receptors, including GPCRs, for which the endogenous ligand has been identified, are referred to as "known" receptors, while receptors for which the endogenous ligand has not been identified are referred to as "orphan" receptors.

GPCRs represent an important area for the development of pharmaceutical products: from approximately 20 of the 100 known GPCRs, approximately 60% of all prescription pharmaceuticals have been developed. For example, in 1999, of the top 100 brand name prescription drugs, the following drugs interact with GPCRs (diseases and/or disorders treated are indicated in parentheses):

15	Claritin® (allergies)	Prozac® (depression)	Vasotec® (hypertension)
	Paxil® (depression)	Zoloft® (depression)	Zyprexa ® (psychotic disorder)
	Cozaar® (hypertension)	Imitrex® (migraine)	Zantac® (reflux)
	Propulsid® (reflux disease)	Risperdal® (schizophrenia)	Serevent® (asthma)
	Pepcid® (reflux)	Gaster® (ulcers)	Atrovent® (bronchospasm)
20	Effexor® (depression)	Depakote® (epilepsy)	Cardura® (prostatic hypertrophy)
	Allegra® (allergies)	Lupron® (prostate cancer)	Zoladex® (prostate cancer)
	Diprivan® (anesthesia)	BuSpar® (anxiety)	Ventolin® (bronchospasm)
	Hytrin® (hypertension)	Wellbutrin® (depression)	Zyrtec® (rhinitis)
	Plavix® (MI/stroke)	Toprol-XL® (hypertension)	Tenormin® (angina)
25	Xalatan® (glaucoma)	Singulair® (asthma)	Diovan® (hypertension)
	Harnal® (prostatic hyperplasia)		

(Med Ad News 1999 Data).

10

. 30

GPCRs share a common structural motif, having seven sequences of between 22 to 24 hydrophobic amino acids that form seven alpha helices, each of which spans the membrane (each span is identified by number, *i.e.*, transmembrane-1 (TM-1), transmebrane-2 (TM-2), *etc.*). The transmembrane helices are joined by strands of amino acids between transmembrane-2 and transmembrane-3, transmembrane-4 and transmembrane-5, and transmembrane-6 and transmembrane-7 on the exterior, or

"extracellular" side, of the cell membrane (these are referred to as "extracellular" regions 1, 2 and 3 (EC-1, EC-2 and EC-3), respectively). The transmembrane helices are also joined by strands of amino acids between transmembrane-1 and transmembrane-2, transmembrane-3 and transmembrane-4, and transmembrane-5 and transmembrane-6 on the interior, or "intracellular" side, of the cell membrane (these are referred to as "intracellular" regions 1, 2 and 3 (IC-1, IC-2 and IC-3), respectively). The "carboxy" ("C") terminus of the receptor lies in the intracellular space within the cell, and the "amino" ("N") terminus of the receptor lies in the extracellular space outside of the cell.

Generally, when an endogenous ligand binds with the receptor (often referred to as "activation" of the receptor), there is a change in the conformation of the intracellular region that allows for coupling between the intracellular region and an intracellular "G-protein." It has been reported that GPCRs are "promiscuous" with respect to G proteins, *i.e.*, that a GPCR can interact with more than one G protein. *See*, Kenakin, T., 43 *Life Sciences* 1095 (1988). Although other G proteins exist, currently, G_q, G_s, G_i, G_z and G_o are G proteins that have been identified. Ligand-activated GPCR coupling with the G-protein initiates a signaling cascade process (referred to as "signal transduction"). Under normal conditions, signal transduction ultimately results in cellular activation or cellular inhibition. Although not wishing to be bound to theory, it is thought that the IC-3 loop as well as the carboxy terminus of the receptor interact with the G protein.

10

20

25

Under physiological conditions, GPCRs exist in the cell membrane in equilibrium between two different conformations: an "inactive" state and an "active" state. A receptor in an inactive state is unable to link to the intracellular signaling transduction pathway to initiate signal transduction leading to a biological response. Changing the receptor conformation to the active state allows linkage to the transduction pathway (via the G-protein) and produces a biological response.

A receptor may be stabilized in an active state by a ligand or a compound such as a drug. Recent discoveries, including but not exclusively limited to modifications to the amino acid sequence of the receptor, provide means other than ligands or drugs to promote and stabilize the receptor in the active state conformation. These means effectively stabilize the receptor in an active state by simulating the effect of a ligand binding to the receptor. Stabilization by such ligand-independent means is termed "constitutive receptor activation."

AREN-0309 SUMMARY OF THE INVENTION

10

15

20

30

PATENT

Disclosed herein are endogenous and non-endogenous versions of human GPCRs and uses thereof.

Some embodiments of the present invention relate to a G protein-coupled receptor encoded by an amino acid sequence of SEQ.ID.NO.:2, non-endogenous, constitutively activated versions of the same, and host cells comprising the same..

Some embodiments of the present invention relate to a plasmid comprising a vector and the cDNA of SEQ.ID.NO.:1 and host cells comprising the same.

Some embodiments of the present invention relate to a G protein-coupled receptor encoded by an amino acid sequence of SEQ.ID.NO.:4, non-endogenous, constitutively activated versions of the same, and host cells comprising the same..

Some embodiments of the present invention relate to a plasmid comprising a vector and the cDNA of SEQ.ID.NO.:3, non-endogenous, constitutively activated versions of the same, and host cells comprising the same.

Some embodiments of the present invention relate to G protein-coupled receptor encoded by an amino acid sequence of SEQ.ID.NO.:6, non-endogenous, constitutively activated versions of the same, and host cells comprising the same.

Some embodiments of the present invention relate to a plasmid comprising a vector and the cDNA of SEQ.ID.NO.:5 and host cells comprising the same.

Some embodiments of the present invention relate to a G protein-coupled receptor encoded by an amino acid sequence of SEQ.ID.NO.:8, non-endogenous, constitutively activated versions of the same, and host cells comprising the same.

Some embodiments of the present invention relate to a plasmid comprising a vector and the cDNA of SEQ.ID.NO.:7, non-endogenous, constitutively activated versions of the same, and host cells comprising the same.

Some embodiments of the present invention relate to a G protein-coupled receptor encoded by an amino acid sequence of SEQ.ID.NO.:10, non-endogenous, constitutively activated versions of the same, and host cells comprising the same.

Some embodiments of the present invention relate to a plasmid comprising a vector and the cDNA of SEQ.ID.NO.:9 and host cells comprising the same.

Some embodiments of the present invention relate to a G protein-coupled receptor encoded by an amino acid sequence of SEQ.ID.NO.:12, non-endogenous, constitutively activated versions of the same, and host cells comprising the same.

Some embodiments of the present invention relate to a plasmid comprising a vector and the cDNA of SEQ.ID.NO.:11, and host cells comprising the same.

Some embodiments of the present invention relate to a G protein-coupled receptor encoded by an amino acid sequence of SEQ.ID.NO.:14, constitutively activated versions of the same, and host cells comprising the same.

Some embodiments of the present invention relate to a plasmid comprising a vector and the cDNA of SEQ.ID.NO.:13 and host cells comprising the same.

Some embodiments of the present invention relate to a G protein-coupled receptor encoded by an amino acid sequence of SEQ.ID.NO.:16, constitutively activated versions of the same, and host cells comprising the same.

Some embodiments of the present invention relate to a plasmid comprising a vector and the cDNA of SEQ.ID.NO.:15 and host cells comprising the same.

Some embodiments of the present invention relate to a G protein-coupled receptor encoded by an amino acid sequence of SEQ.ID.NO.:18, constitutively activated versions of the same, and host cells comprising the same.

Some embodiments of the present invention relate to a plasmid comprising a vector and the cDNA of SEQ.ID.NO.:17 and host cells comprising the same.

Some embodiments of the present invention relate to a G protein-coupled receptor encoded by an amino acid sequence of SEQ.ID.NO.:20, constitutively activated versions of the same, and host cells comprising the same.

Some embodiments of the present invention relate to a plasmid comprising a vector and the cDNA of SEQ.ID.NO.:19 and host cells comprising the same.

BRIEF DESCRIPTION OF THE DRAWINGS

15

20

30

Figure 1 is a graphic representation of activation of RUP32, $G_q(del)/G_i$, RUP32 co-transfected with $G_q(del)/G_i$, and CMV (control; expression vector) in a second messenger assay measuring the accumulation of inositol phosphate (IP₃) utilizing 293 cells.

Figure 2 provides an illustration of second messenger IP₃ production from endogenous version RUP35 and RUP36 as compared with the control ("CMV").

DETAILED DESCRIPTION

10

The scientific literature that has evolved around receptors has adopted a number of terms to refer to ligands having various effects on receptors. For clarity and consistency, the following definitions will be used throughout this patent document. To the extent that these definitions conflict with other definitions for these terms, the following definitions shall control:

AGONISTS shall mean materials (e.g., ligands, candidate compounds) that activate the intracellular response when they bind to the receptor, or enhance GTP binding to membranes. In some embodiments, AGONISTS are those materials not previously known to activate the intracellular response when they bind to the receptor or to enhance GTP binding to membranes.

AMINO ACID ABBREVIATIONS used herein are set out in Table A:

TABLE A

ALANINE	ALA	Α .
ARGININE	ARG	R
ASPARAGINE	ASN	N
ASPARTIC ACID	ASP	. D
CYSTEINE	CYS	С
GLUTAMIC ACID	GLU	Е
GLUTAMINE	GLN	Q
GLYCINE	GLY	G
HISTIDINE	HIS	Н
ISOLEUCINE	ILE	I
LEUCINE	LEU	L
LYSINE	LYS	K

AREN-0309 PATENT

MET	М
PHE	F
PRO	Р .
SER	S
THR	T
TRP	W
TYR	Y
VAL	V
	PHE PRO SER THR TRP TYR

ANTAGONIST shall mean materials (e.g., ligands, candidate compounds) that competitively bind to the receptor at the same site as the agonists but which do not activate the intracellular response initiated by the active form of the receptor, and can thereby inhibit the intracellular responses by agonists. ANTAGONISTS do not diminish the baseline intracellular response in the absence of an agonist. In some embodiments, ANTAGONISTS are those materials not previously known to activate the intracellular response when they bind to the receptor or to enhance GTP binding to membranes.

CANDIDATE COMPOUND shall mean a molecule (for example, and not limitation, a chemical compound) that is amenable to a screening technique. Preferably, the phrase "candidate compound" does not include compounds which were publicly known to be compounds selected from the group consisting of inverse agonist, agonist or antagonist to a receptor, as previously determined by an indirect identification process ("indirectly identified compound"); more preferably, not including an indirectly identified compound which has previously been determined to have therapeutic efficacy in at least one mammal; and, most preferably, not including an indirectly identified compound which has previously been determined to have therapeutic utility in humans.

COMPOSITION means a material comprising at least one component; a "pharmaceutical composition" is an example of a composition.

COMPOUND EFFICACY shall mean a measurement of the ability of a compound to inhibit or stimulate receptor functionality; i.e. the ability to activate/inhibit a

20

signal transduction pathway, as opposed to receptor binding affinity. Exemplary means of detecting compound efficacy are disclosed in the Example section of this patent document.

CODON shall mean a grouping of three nucleotides (or equivalents to nucleotides)

which generally comprise a nucleoside (adenosine (A), guanosine (G), cytidine (C), uridine

(U) and thymidine (T)) coupled to a phosphate group and which, when translated, encodes an amino acid.

CONSTITUTIVELY ACTIVATED RECEPTOR shall mean a receptor subjected to constitutive receptor activation. A constitutively activated receptor can be endogenous or non-endogenous.

10

15

25

CONSTITUTIVE RECEPTOR ACTIVATION shall mean stabilization of a receptor in the active state by means other than binding of the receptor with its ligand or a chemical equivalent thereof.

CONTACT or CONTACTING shall mean bringing at least two moieties together, whether in an in vitro system or an in vivo system.

DIRECTLY IDENTIFYING or DIRECTLY IDENTIFIED, in relationship to the phrase "candidate compound", shall mean the screening of a candidate compound against a constitutively activated receptor, preferably a constitutively activated orphan receptor, and most preferably against a constitutively activated G protein-coupled cell surface orphan receptor, and assessing the compound efficacy of such compound. This phrase is, under no circumstances, to be interpreted or understood to be encompassed by or to encompass the phrase "indirectly identifying" or "indirectly identified."

ENDOGENOUS shall mean a material that a mammal naturally produces. ENDOGENOUS in reference to, for example and not limitation, the term "receptor," shall mean that which is naturally produced by a mammal (for example, and not limitation, a human) or a virus. By contrast, the term NON-ENDOGENOUS in this context shall mean that which is not naturally produced by a mammal (for example, and not limitation, a human) or a virus. For example, and not limitation, a receptor which is not constitutively active in its endogenous form, but when manipulated becomes constitutively active, is most preferably referred to herein as a "non-endogenous, constitutively activated receptor." Both terms can be utilized to describe both "in vivo" and "in vitro" systems. For example, and not limitation, in a screening approach, the endogenous or non-endogenous receptor may be

5

20

25

in reference to an in vitro screening system. As a further example and not limitation, where the genome of a mammal has been manipulated to include a non-endogenous constitutively activated receptor, screening of a candidate compound by means of an in vivo system is viable.

G PROTEIN COUPLED RECEPTOR FUSION PROTEIN and **GPCR FUSION PROTEIN**, in the context of the invention disclosed herein, each mean a non-endogenous protein comprising an endogenous, constitutively activate GPCR or a non-endogenous, constitutively activated GPCR fused to at least one G protein, most preferably the alpha (α) subunit of such G protein (this being the subunit that binds GTP), with the G protein preferably being of the same type as the G protein that naturally couples with endogenous orphan GPCR. For example, and not limitation, in an endogenous state, if the G protein " $G_s\alpha$ " is the predominate G protein that couples with the GPCR, a GPCR Fusion Protein based upon the specific GPCR would be a non-endogenous protein comprising the GPCR fused to $G_s\alpha$; in some circumstances, as will be set forth below, a non-predominant G protein can be fused to the GPCR or there may be spacers between the two.

HOST CELL shall mean a cell capable of having a Plasmid and/or Vector incorporated therein. In the case of a prokaryotic Host Cell, a Plasmid is typically replicated as a autonomous molecule as the Host Cell replicates (generally, the Plasmid is thereafter isolated for introduction into a eukaryotic Host Cell); in the case of a eukaryotic Host Cell, a Plasmid is integrated into the cellular DNA of the Host Cell such that when the eukaryotic Host Cell replicates, the Plasmid replicates. In some embodiments the Host Cell is eukaryotic, more preferably, mammalian, and most preferably selected from the group consisting of 293, 293T and COS-7 cells.

INDIRECTLY IDENTIFYING or INDIRECTLY IDENTIFIED means the traditional approach to the drug discovery process involving identification of an endogenous ligand specific for an endogenous receptor, screening of candidate compounds against the receptor for determination of those which interfere and/or compete with the ligand-receptor interaction, and assessing the efficacy of the compound for affecting at least one second messenger pathway associated with the activated receptor.

INHIBIT or **INHIBITING**, in relationship to the term "response" shall mean that a response is decreased or prevented in the presence of a compound as opposed to in the absence of the compound.

INVERSE AGONISTS shall mean materials (e.g., ligand, candidate compound) which bind to either the endogenous form of the receptor or to the constitutively activated form of the receptor, and which inhibit the baseline intracellular response initiated by the active form of the receptor below the normal base level of activity which is observed in the absence of agonists, or decrease GTP binding to membranes. Preferably, the baseline intracellular response is inhibited in the presence of the inverse agonist by at least 30%, at least 50%, at least 60%, at least 70%, at least 75%, at least 80%, at least 85%, at least 90%, at least 92%, at least 94%, at least 95%, at least 97%, at least 98%, and most preferably at least 99% as compared with the baseline response in the absence of the inverse agonist.

KNOWN RECEPTOR shall mean an endogenous receptor for which the endogenous ligand specific for that receptor has been identified.

15

20

25

30

LIGAND shall mean a molecule specific for a naturally occurring receptor.

MUTANT or MUTATION in reference to an endogenous receptor's nucleic acid and/or amino acid sequence shall mean a specified change or changes to such endogenous sequences such that a mutated form of an endogenous, non-constitutively activated receptor evidences constitutive activation of the receptor. In terms of equivalents to specific sequences, a subsequent mutated form of a human receptor is considered to be equivalent to a first mutation of the human receptor if (a) the level of constitutive activation of the subsequent mutated form of a human receptor is substantially the same as that evidenced by the first mutation of the receptor; and (b) the percent sequence (amino acid and/or nucleic acid) homology between the subsequent mutated form of the receptor and the first mutation of the receptor is at least 80%, at least 85%, at least 90%, at least 92%, at least 94%, at least 95%, at least 96%, at least 97%, at least 98%, and most preferably at least 99%. In some embodiments, owing to the fact that some preferred cassettes disclosed herein for achieving constitutive activation include a single amino acid and/or codon change between the endogenous and the non-endogenous forms of the GPCR, it is preferred that the percent sequence homology should be at least 98%.

NON-ORPHAN RECEPTOR shall mean an endogenous naturally occurring molecule specific for an identified ligand wherein the binding of a ligand to a receptor activates an intracellular signaling pathway.

ORPHAN RECEPTOR shall mean an endogenous receptor for which the ligand specific for that receptor has not been identified or is not known.

PHARMACEUTICAL COMPOSITION shall mean a composition comprising at least one active ingredient, whereby the composition is amenable to investigation for a specified, efficacious outcome in a mammal (for example, and not limitation, a human). Those of ordinary skill in the art will understand and appreciate the techniques appropriate for determining whether an active ingredient has a desired efficacious outcome based upon the needs of the artisan.

PLASMID shall mean the combination of a Vector and cDNA. Generally, a Plasmid is introduced into a Host Cell for the purposes of replication and/or expression of the cDNA as a protein.

15

20

25

SECOND MESSENGER shall mean an intracellular response produced as a result of receptor activation. A second messenger can include, for example, inositol triphosphate (IP₃), diacycglycerol (DAG), cyclic AMP (cAMP), and cyclic GMP (cGMP). Second messenger response can be measured for a determination of receptor activation. In addition, second messenger response can be measured for the direct identification of candidate compounds, including for example, inverse agonists, agonists, and antagonists.

SIGNAL TO NOISE RATIO shall mean the signal generated in response to activation, amplification, or stimulation wherein the signal is above the background noise or the basal level in response to non-activation, non-amplification, or non-stimulation.

SPACER shall mean a translated number of amino acids that are located after the last codon or last amino acid of a gene, for example a GPCR of interest, but before the start codon or beginning regions of the G protein of interest, wherein the translated number amino acids are placed in-frame with the beginnings regions of the G protein of interest. The number of translated amino acids can be tailored according to the needs of the skilled artisan and is generally from about one amino acid, preferably two amino acids, more preferably three amino acids, more preferably four amino acids, more preferably six amino acids, more preferably seven amino acids, more preferably

eight amino acids, more preferably nine amino acids, more preferably ten amino acids, more preferably eleven amino acids, and even more preferably twelve amino acids.

STIMULATE or STIMULATING, in relationship to the term "response" shall mean that a response is increased in the presence of a compound as opposed to in the absence of the compound.

SUBSTANTIALLY shall refer to a result which is within 40% of a control result, preferably within 35%, more preferably within 30%, more preferably within 25%, more preferably within 20%, more preferably within 15%, more preferably within 10%, more preferably within 5%, more preferably within 2%, and most preferably within 1% of a control result. For example, in the context of receptor functionality, a test receptor may exhibit substantially similar results to a control receptor if the transduced signal, measured using a method taught herein or similar method known to the art-skilled, if within 40% of the signal produced by a control signal.

VECTOR in reference to cDNA shall mean a circular DNA capable of incorporating at least one cDNA and capable of incorporation into a Host Cell.

The order of the following sections is set forth for presentational efficiency and is not intended, nor should be construed, as a limitation on the disclosure or the claims to follow.

20 A. Introduction

15

30

The traditional study of receptors has typically proceeded from the a priori assumption (historically based) that the endogenous ligand must first be identified before discovery could proceed to find antagonists and other molecules that could affect the receptor. Even in cases where an antagonist might have been known first, the search 25 immediately extended to looking for the endogenous ligand. This mode of thinking has persisted in receptor research even after the discovery of constitutively activated receptors. What has not been heretofore recognized is that it is the active state of the receptor that is most useful for discovering agonists and inverse agonists of the receptor. For those diseases which result from an overly active receptor or an under-active receptor, what is desired in a therapeutic drug is a compound which acts to diminish the active state of a receptor or enhance the activity of the receptor, respectively, not necessarily a drug which is an antagonist to the endogenous ligand. This is because a compound that reduces or enhances

the activity of the active receptor state need not bind at the same site as the endogenous ligand. Thus, as taught by a method of this invention, any search for therapeutic compounds should start by screening compounds against the ligand-independent active state.

5 B. Identification of Human GPCRs

10

15

The efforts of the Human Genome project has led to the identification of a plethora of information regarding nucleic acid sequences located within the human genome; it has been the case in this endeavor that genetic sequence information has been made available without an understanding or recognition as to whether or not any particular genomic sequence does or may contain open-reading frame information that translate human proteins. Several methods of identifying nucleic acid sequences within the human genome are within the purview of those having ordinary skill in the art. For example, and not limitation, a variety of human GPCRs, disclosed herein, were discovered by reviewing the GenBankTM database. Table B, below, lists several endogenous GPCRs that we have discovered, along with other GPCRs that are homologous to the disclosed GPCR.

TABLE B

Disclosed	Accession	Open Reading	Reference To	Per Cent
Human	Number	Frame	Homologous	Homology
Orphan GPCRs	Identified	(Base Pairs)	GPCR	To Designated
				GPCR
hRUP28	AC073957	1,002bp	hGPR30	34%
hRUP29	AC083865	918bp	hGPR18	27%
hRUP30	AC055863	1,125bp	hBRB1	27%
hRUP31	AL356214	1,086bp	hGALR-1	31%
hRUP32	AL513524	1,038bp	hPNR	43%
hRUP33	AL513524	1,020bp	GPR57	50%
			GPR58	51%
hRUP34	AL513524	1,029bp	hPNR	45%
hRUP35	AC021089	1,062bp	hк-type 3	27%
			opioid	,
hRUP36	AC090099	969bp	GPR90	42%

AREN-0309				PATENT	
hRUP37	AC090099	969bp	hMRG	41%	

Receptor homology is useful in terms of gaining an appreciation of a role of the receptors within the human body. As the patent document progresses, techniques for mutating these receptors to establish non-endogenous, constitutively activated versions of these receptors will be discussed.

The techniques disclosed herein are also applicable to other human GPCRs known to the art, as will be apparent to those skilled in the art.

C. Receptor Screening

15

20

25

30

. ____

Screening candidate compounds against a non-endogenous, constitutively activated version of the GPCRs disclosed herein allows for the direct identification of candidate compounds which act at the cell surface receptor, without requiring use of the receptor's endogenous ligand. Using routine, and often commercially available techniques, one can determine areas within the body where the endogenous version of human GPCRs disclosed herein is expressed and/or over-expressed. The expression location of a receptor in a specific tissue provides a scientist with the ability to assign a physiological functional role of the receptor. It is also possible using these techniques to determine related disease/disorder states which are associated with the expression and/or over-expression of the receptor; such an approach is disclosed in this patent document. Furthermore, expression of a receptor in diseased organs can assist one in determining the magnitude of the clinical relevance of the receptor.

Constitutive activation of the GPCRs disclosed herein is based upon the distance from the proline residue at which is presumed to be located within TM6 of the GPCR; this algorithmic technique is disclosed in co-pending and commonly assigned patent document PCT Application Number PCT/US99/23938, published as WO 00/22129 on April 20, 2000, which, along with the other patent documents listed herein, is incorporated herein by reference. The algorithmic technique is not predicated upon traditional sequence "alignment" but rather a specified distance from the aforementioned TM6 proline residue (or, of course, endogenous constitutive substitution for such proline residue). By mutating the amino acid residue located 16 amino acid residues from this residue (presumably located in the IC3 region of the receptor) to, most preferably, a lysine residue, constitutive

activation of the receptor may be obtained. Other amino acid residues may be useful in the mutation at this position to achieve this objective and will be discussed in detail, below.

D. Disease/Disorder Identification and/or Selection

20

25

As will be set forth in greater detail below, inverse agonists and agonists to the non-endogenous, constitutively activated GPCR can be identified by the methodologies of this invention. Such inverse agonists and agonists are ideal candidates as lead compounds in drug discovery programs for treating diseases related to this receptor. Because of the ability to directly identify inverse agonists to the GPCR, thereby allowing for the development of pharmaceutical compositions, a search for diseases and disorders associated with the GPCR is relevant. The expression location of a receptor in a specific tissue provides a scientist with the ability to assign a physiological function to the receptor. For example, scanning both diseased and normal tissue samples for the presence of the GPCR now becomes more than an academic exercise or one which might be pursued along the path of identifying an endogenous ligand to the specific GPCR. Tissue scans can be conducted across a broad range of healthy and diseased tissues. Such tissue scans provide a potential first step in associating a specific receptor with a disease and/or disorder. Furthermore, expression of a receptor in diseased organs can assist one in determining the magnitude of clinical relevance of the receptor.

The DNA sequence of the GPCR can be used to make a probe/primer. In some preferred embodiments the DNA sequence is used to make a probe for (a) dot-blot analysis against tissue-mRNA, and/or (b) RT-PCR identification of the expression of the receptor in tissue samples. The presence of a receptor in a tissue source, or a diseased tissue, or the presence of the receptor at elevated concentrations in diseased tissue compared to a normal tissue, can be used to correlate location to function and indicate the receptor's physiological role/function and create a treatment regimen, including but not limited to, a disease associated with that function/role. Receptors can also be localized to regions of organs by this technique. Based on the known or assumed roles/functions of the specific tissues to which the receptor is localized, the putative physiological function of the receptor can be deduced. For example and not limitation, proteins located/expressed in areas of the thalamus are associated with sensorimotor processing and arousal (see, Goodman & Gilman's, The Pharmacological Basis of Therapeutics, 9th Edition, page 465 (1996)). Proteins expressed in the hippocampus or in Schwann cells are associated with learning and

memory, and myelination of peripheral nerves, respectively (*see*, Kandel, E. et al., Essentials of Neural Science and Behavior pages 657, 680 and 28, respectively (1995)).

E. Screening of Candidate Compounds

5

20

25

1. Generic GPCR screening assay techniques

When a G protein receptor becomes constitutively active, it binds to a G protein (e.g., G_q, G_s, G_i, G_z, Go) and stimulates the binding of GTP to the G protein. The G protein then acts as a GTPase and hydrolyzes the GTP to GDP, whereby the receptor, under normal conditions, becomes deactivated. However, constitutively activated receptors continue to exchange GDP to GTP. A non-hydrolyzable analog of GTP, [35S]GTPγS, can be used to monitor enhanced binding to membranes which express constitutively activated receptors. It is reported that [35S]GTPγS can be used to monitor G protein coupling to membranes in the absence and presence of ligand. An example of this monitoring, among other examples well-known and available to those in the art, was reported by Traynor and Nahorski in 1995. The use of this assay system is typically for initial screening of candidate compounds because the system is generically applicable to all G protein-coupled receptors regardless of the particular G protein that interacts with the intracellular domain of the receptor.

2. Specific GPCR screening assay techniques

Once candidate compounds are identified using the "generic" G protein-coupled receptor assay (i.e., an assay to select compounds that are agonists or inverse agonists), further screening to confirm that the compounds have interacted at the receptor site is preferred. For example, a compound identified by the "generic" assay may not bind to the receptor, but may instead merely "uncouple" the G protein from the intracellular domain.

a. G_{s} , G_{z} and G_{i} .

G_s stimulates the enzyme adenylyl cyclase. G_i (and G_z and Go), on the other hand, inhibits adenylyl cyclase. Adenylyl cyclase catalyzes the conversion of ATP to cAMP; thus, constitutively activated GPCRs that couple the G_s protein are associated with increased cellular levels of cAMP. On the other hand, constitutively activated GPCRs that couple G_i (or G_z, Go) protein are associated with decreased cellular levels of cAMP. See, generally, "Indirect Mechanisms of Synaptic Transmission," Chpt. 8, From Neuron To Brain (3rd Ed.) Nichols, J.G. et al eds. Sinauer Associates, Inc. (1992). Thus, assays that detect cAMP can be utilized to determine if a candidate compound is, e.g., an inverse

agonist to the receptor (*i.e.*, such a compound would decrease the levels of cAMP). A variety of approaches known in the art for measuring cAMP can be utilized; a most preferred approach relies upon the use of anti-cAMP antibodies in an ELISA-based format. Another type of assay that can be utilized is a whole cell second messenger reporter system assay. Promoters on genes drive the expression of the proteins that a particular gene encodes. Cyclic AMP drives gene expression by promoting the binding of a cAMP-responsive DNA binding protein or transcription factor (CREB) that then binds to the promoter at specific sites (cAMP response elements) and drives the expression of the gene. Reporter systems can be constructed which have a promoter containing multiple cAMP response elements before the reporter gene, *e.g.*, β-galactosidase or luciferase. Thus, a constitutively activated G_s-linked receptor causes the accumulation of cAMP that then activates the gene and leads to the expression of the reporter protein. The reporter protein such as β-galactosidase or luciferase can then be detected using standard biochemical assays (Chen et al. 1995).

b. G_o and G_{qr}

15

25

30

G_q and G_o are associated with activation of the enzyme phospholipase C, which in turn hydrolyzes the phospholipid PIP₂, releasing two intracellular messengers: diacycloglycerol (DAG) and inositol 1,4,5-triphoisphate (IP₃). Increased accumulation of IP₃ is associated with activation of G_q- and Go-associated receptors. See, generally, "Indirect Mechanisms of Synaptic Transmission," Chpt. 8, From Neuron To Brain (3rd Ed.) Nichols, J.G. et al eds. Sinauer Associates, Inc. (1992). Assays that detect IP₃ accumulation can be utilized to determine if a candidate compound is, e.g., an inverse agonist to a G_q- or Go-associated receptor (i.e., such a compound would decrease the levels of IP₃). G_q- associated receptors can also be examined using an AP1 reporter assay wherein G_q-dependent phospholipase C causes activation of genes containing AP1 elements; thus, activated G_q-associated receptors will evidence an increase in the expression of such genes, whereby inverse agonists thereto will evidence a decrease in such expression, and agonists will evidence an increase in such expression. Commercially available assays for such detection are available.

3. GPCR Fusion Protein

The use of an endogenous, constitutively activated GPCR or a non-endogenous, constitutively activated GPCR, for use in screening of candidate compounds for the direct identification of inverse agonists, agonists provide an interesting screening challenge in that, by definition, the receptor is active even in the absence of an endogenous ligand bound thereto. Thus, in order to differentiate between, e.g., the non-endogenous receptor in the presence of a candidate compound and the non-endogenous receptor in the absence of that compound, with an aim of such a differentiation to allow for an understanding as to whether such compound may be an inverse agonist or agonist or have no affect on such a receptor, it is preferred that an approach be utilized that can enhance such differentiation. A preferred approach is the use of a GPCR Fusion Protein.

10

20

25

30

Generally, once it is determined that a non-endogenous GPCR has been constitutively activated using the assay techniques set forth above (as well as others), it is possible to determine the predominant G protein that couples with the endogenous GPCR. Coupling of the G protein to the GPCR provides a signaling pathway that can be assessed. In some embodiments it is preferred that screening take place using a mammalian expression system, such a system will be expected to have endogenous G protein therein. Thus, by definition, in such a system, the non-endogenous, constitutively activated GPCR will continuously signal. In some embodiments it is preferred that this signal be enhanced such that in the presence of, e.g., an inverse agonist to the receptor, it is more likely that it will be able to more readily differentiate, particularly in the context of screening, between the receptor when it is contacted with the inverse agonist.

The GPCR Fusion Protein is intended to enhance the efficacy of G protein coupling with the non-endogenous GPCR. The GPCR Fusion Protein is preferred for screening with either an endogenous, constitutively active GPCR or a non-endogenous, constitutively activated GPCR because such an approach increases the signal that is utilized in such screening techniques. This is important in facilitating a significant "signal to noise" ratio; such a significant ratio is preferred for the screening of candidate compounds as disclosed herein.

The construction of a construct useful for expression of a GPCR Fusion Protein is within the purview of those having ordinary skill in the art. Commercially available expression vectors and systems offer a variety of approaches that can fit the particular needs of an investigator. Important criteria on the construction of such a GPCR Fusion Protein

construct include but are not limited to, that the endogenous GPCR sequence and the G protein sequence both be in-frame (preferably, the sequence for the endogenous GPCR is upstream of the G protein sequence), and that the "stop" codon of the GPCR be deleted or replaced such that upon expression of the GPCR, the G protein can also be expressed. Other embodiments include constructs wherein the endogenous GPCR sequence and the G protein sequence are not in-frame and/or the "stop" codon is not deleted or replaced. The GPCR can be linked directly to the G protein, or there can be spacer residues between the two (preferably, no more than about 12, although this number can be readily ascertained by one of ordinary skill in the art). Based upon convenience it is preferred to use a spacer. Preferably, the G protein that couples to the non-endogenous GPCR will have been identified prior to the creation of the GPCR Fusion Protein construct. Because there are only a few G proteins that have been identified, it is preferred that a construct comprising the sequence of the G protein (i.e., a universal G protein construct (see Examples)) be available for insertion of an endogenous GPCR sequence therein; this provides for further efficiency in the context of large-scale screening of a variety of different endogenous GPCRs having different sequences.

15

20

30

As noted above, constitutively activated GPCRs that couple to G_i , G_z and G_o are expected to inhibit the formation of cAMP making assays based upon these types of GPCRs challenging (*i.e.*, the cAMP signal decreases upon activation thus making the direct identification of, *e.g.*, inverse agonists (which would further decrease this signal), challenging. As will be disclosed herein, we have ascertained that for these types of receptors, it is possible to create a GPCR Fusion Protein that is not based upon the GPCRs endogenous G protein, in an effort to establish a viable cyclase-based assay. Thus, for example, an endogenous G_i coupled receptor can be fused to a G_s protein—such a fusion construct, upon expression, "drives" or "forces" the endogenous GPCR to couple with, *e.g.*, G_s rather than the "natural" G_i protein, such that a cyclase-based assay can be established. Thus, for G_i , G_z and G_o coupled receptors, in some embodiments it is preferred that when a GPCR Fusion Protein is used and the assay is based upon detection of adenylyl cyclase activity, that the fusion construct be established with G_s (or an equivalent G protein that stimulates the formation of the enzyme adenylyl cyclase).

	, 			
1 C	Effect of cAMP	Effect of IP ₃	Tresont of	Effect on IP ₃
· ·	Effect of CAIVE	Lifect of IL'3	Effect of	Enect on 112
1		•		

protein	Production upon Activation of GPCR (i.e., constitutive activation or agonist binding)	Accumulation upon Activation of GPCR (i.e., constitutive activation or agonist binding)	cAMP Production upon contact with an Inverse Agonist	Accumulation upon contact with an Inverse Agonist
Gs	Increase	N/A	Decrease	N/A
Gi	Decrease	N/A	Increase	N/A
Gz	Decrease	N/A	Increase	N/A
Go	Decrease	Increase	Increase	Decrease
G_q	N/A	Increase	N/A	Decrease

Equally effective is a G Protein Fusion construct that utilizes a G_q Protein fused with a G_s , G_i , G_z or G_o Protein. In some embodiments a preferred fusion construct can be accomplished with a G_q Protein wherein the first six (6) amino acids of the G-protein α -subunit (" $G\alpha q$ ") is deleted and the last five (5) amino acids at the C-terminal end of $G\alpha q$ is replaced with the corresponding amino acids of the $G\alpha$ of the G protein of interest. For example, a fusion construct can have a G_q (6 amino acid deletion) fused with a G_i Protein, resulting in a " G_q/G_i Fusion Construct". This fusion construct will forces the endogenous G_i coupled receptor to couple to its non-endogenous G protein, G_q , such that the second messenger, for example, inositol triphosphate or diacylgycerol, can be measured *in lieu* of cAMP production.

4. Co-transfection of a Target G_i Coupled GPCR with a Signal-Enhancer G_s Coupled GPCR (cAMP Based Assays)

15

A G_i coupled receptor is known to inhibit adenylyl cyclase, and, therefore, decreases the level of cAMP production, which can make assessment of cAMP levels challenging. An effective technique in measuring the decrease in production of cAMP as an indication of constitutive activation of a receptor that predominantly couples G_i upon activation can be accomplished by co-transfecting a signal enhancer, e.g., a non-endogenous, constitutively activated receptor that predominantly couples with G_s upon activation (e.g., TSHR-A623I, disclosed below), with the G_i linked GPCR. As is apparent, constitutive activation of a G_s

coupled receptor can be determined based upon an increase in production of cAMP. Constitutive activation of a G_i coupled receptor leads to a decrease in production cAMP. Thus, the co-transfection approach is intended to advantageously exploit these "opposite" affects. For example, co-transfection of a non-endogenous, constitutively activated G_s coupled receptor (the "signal enhancer") with the endogenous G_i coupled receptor (the "target receptor") provides a baseline cAMP signal (*i.e.*, although the G_i coupled receptor will decrease cAMP levels, this "decrease" will be relative to the substantial increase in cAMP levels established by constitutively activated G_s coupled signal enhancer). By then co-transfecting the signal enhancer with a constitutively activated version of the target receptor, cAMP would be expected to further decrease (relative to base line) due to the increased functional activity of the G_i target (*i.e.*, which decreases cAMP).

Screening of candidate compounds using a cAMP based assay can then be accomplished, with two 'changes' relative to the use of the endogenous receptor/G-protein fusion: first, relative to the G_i coupled target receptor, "opposite" effects will result, *i.e.*, an inverse agonist of the G_i coupled target receptor will increase the measured cAMP signal, while an agonist of the G_i coupled target receptor will decrease this signal; second, as would be apparent, candidate compounds that are directly identified using this approach should be assessed independently to ensure that these do not target the signal enhancing receptor (this can be done prior to or after screening against the co-transfected receptors).

20

25

15

F. Medicinal Chemistry

Generally, but not always, direct identification of candidate compounds is conducted in conjunction with compounds generated via combinatorial chemistry techniques, whereby thousands of compounds are randomly prepared for such analysis. Generally, the results of such screening will be compounds having unique core structures; thereafter, these compounds may be subjected to additional chemical modification around a preferred core structure(s) to further enhance the medicinal properties thereof. Such techniques are known to those in the art and will not be addressed in detail in this patent document.

30

G. Pharmaceutical compositions

Candidate compounds selected for further development can be formulated into pharmaceutical compositions using techniques well known to those in the art. Suitable pharmaceutically-acceptable carriers are available to those in the art; for example, see Remington's Pharmaceutical Sciences, 16th Edition, 1980, Mack Publishing Co., (Osol et al., eds.).

H. Other Utilities

10

20

25

Although a preferred use of the non-endogenous versions of the GPCRs disclosed herein may be for the direct identification of candidate compounds as inverse agonists or agonists (preferably for use as pharmaceutical agents), other uses of these versions of GPCRs exist. For example, in vitro and in vivo systems incorporating GPCRs can be utilized to further elucidate and understand the roles these receptors play in the human condition, both normal and diseased, as well as understanding the role of constitutive activation as it applies to understanding the signaling cascade. In some embodiments it is preferred that the endogenous receptors be "orphan receptors", i.e., the endogenous ligand for the receptor has not been identified. In some embodiments, therefore, the modified, non-endogenous GPCRs can be used to understand the role of endogenous receptors in the human body before the endogenous ligand therefore is identified. Such receptors can also be used to further elucidate known receptors and the pathways through which they transduce a signal. Other uses of the disclosed receptors will become apparent to those in the art based upon, inter alia, a review of this patent document.

EXAMPLES

The following examples are presented for purposes of elucidation, and not limitation, of the present invention. While specific nucleic acid and amino acid sequences are disclosed herein, those of ordinary skill in the art are credited with the ability to make minor modifications to these sequences while achieving the same or substantially similar results reported below. The traditional approach to application or understanding of sequence cassettes from one sequence to another (e.g. from rat receptor to human receptor or from human receptor A to human receptor B) is generally predicated upon sequence alignment techniques whereby the sequences are aligned in an effort to determine areas of commonality. The mutational approach disclosed herein does not rely upon this approach but is instead based upon an algorithmic approach and a positional distance from a conserved proline residue located within the TM6 region of human GPCRs. Once this

approach is secured, those in the art are credited with the ability to make minor modifications thereto to achieve substantially the same results (i.e., constitutive activation) disclosed herein. Such modified approaches are considered within the purview of this disclosure.

5 Example 1

ENDOGENOUS HUMAN GPCRS

1. Identification of Human GPCRs

The disclosed endogenous human GPCRs were identified based upon a review of the GenBankTM database information. While searching the database, the following cDNA clones were identified as evidenced below (Table C).

Disclosed Accession Open Reference To Nucleic Acid Amino Acid Human Number Reading Homologous SEQ.ID. NO. SEQ.ID.NO. Orphan Identified Frame **GPCR GPCRs** (Base Pairs) hRUP28 AC073957 1,002bp hGPR30 ī 2 hRUP29 AC083865 918bp hGPR18 3 4 hRUP30 AC055863 1,125bp hBRB1 5 6 hRUP31 AL356214 1,086bp 7 8 hGALR-1 hRUP32 1,038bp AL513524 hPNR 9 10 hRUP33 AL513524 1,020bp GPR57 11 12 GPR58 hRUP34 AL513524 1,029bp 13 hPNR 14 hRUP35 AC021089 1,062bp hк-type 3 15 16 opioid hRUP36 AC090099 969bp GPR90 17 18 hRUP37 AC090099 969bp hMRG 19 20

TABLE C

2. Full Length Cloning

a. hRUP28 (Seq. Id. Nos. 1 & 2)

15 The disclosed human RUP28 was identified based upon the use of GenBank database information. While searching the database, a cDNA clone with Accession Number AC073957 was identified as a human genomic sequence from chromosome 7.

The full length RUP28 was cloned by PCR using primers:

15

- 3' (SEQ.ID.NO.:21; sense, 5' of initiation codon),
- 5'-CTGCGTCCACCAGAGTCACGTCTCC-3' (SEQ.ID.NO.:22; antisense, 3' of stop codon), and human adult liver Marathon-ReadyTM cDNA (Clontech) as template.
- 5 AdvantageTM cDNA polymerase (Clontech) was used for the amplification in a 50μl reaction by the following cycle with step 2 to 4 repeated 35 times: 95°C for 5 min; 94°C for 30 sec; 58°C for 30 sec; 72°C for 1 min 30 sec; and 72°C for 7 min.

A 1.16kb PCR fragment was isolated from a 1% agarose gel and cloned into the pCRII-TOPO vector (Invitrogen) and sequenced using the ABI Big Dye Terminator Kit (P.E. Biosystems). See, SEQ.ID.NO.:1 for the nucleic acid sequence and SEQ.ID.NO.:2 for the putative amino acid sequence.

b. hRUP29 (Seq. Id. Nos. 3 & 4)

The disclosed human RUP29 was identified based upon the use of GenBank database information. While searching the database, a cDNA clone with Accession Number AC0083865 was identified as a human genomic sequence from chromosome 7.

The full length RUP29 was cloned by PCR using primers:

- 5'-GTATGCCTGGCCACAATACCTCCAGG-3' (SEQ.ID.NO.:23; sense, containing the initiation codon),
- 5'-GTTTGTGGCTAACGGCACAAAACACAATTCC-3' (SEQ.ID.NO.:24; antisense, containing the stop codon) and human genomic DNA as template. TaqPlus[®] Precision DNA polymerase (Stratagene) was used for the amplification in a 50µl reaction by the following cycle with step 2 to 4 repeated 35 times: 94°C for 5 min; 94°C for 30 sec; 54°C for 30 sec; 72°C for 1 min 30 sec; and 72°C for 7 min.

A 930bp PCR fragment was isolated from a 1% agarose gel and cloned into the pCRII-TOPO vector (Invitrogen) and sequenced using the ABI Big Dye Terminator Kit (P.E. Biosystems).

Rapid amplification of cDNA ends (RACE) was performed using human leukocyte and ovary Marathon-ReadyTM cDNA (Clontech) to determine the precise 5' end of RUP29 cDNA. RUP29 specific primer (1) having the sequence:

30 5-GGTACCACAATGACAATCACCAGCGTCC-3'(SEQ.ID.NO.:25) and AP1 primer (Clontech) were used for the first-round PCR reaction, and RUP29 specific primer (2) having the following sequence:

5'-

15

GGAACGTGAGGTACATGTGGATGTGCAGC-3' (SEQ.ID.NO.:26)

and AP2 primer (Clontech) were used for the second-round PCR reaction. The products of the RACE reactions were isolated and cloned into the pCRII-TOPO vector (Invitrogen) and sequenced. See, SEQ.ID.NO.:3 for the nucleic acid sequence and SEQ.ID.NO.:4 for the putative amino acid sequence.

c. hRUP30 (Seq. Id. Nos. 5 & 6)

The disclosed human RUP30 was identified based upon the use of GenBank database information. While searching the database, a cDNA clone with Accession Number AC055863 was identified as a human genomic sequence from chromosome 17.

The full length RUP30 was cloned by 5'RACE -PCR with a human pancreas Marathon-ReadyTM cDNA (Clontech) as template and the following oligonucleotide:

- 5'-GCAGTGTAGCGGTCAACCGTGAGCAGG-3'(SEQ.ID.NO.:27; sense, containing the initiation codon), and AP1 primer (Clontech) were used for the first round of RT-PCR and oligonucleotide:
- 5'-TGAGCAGGATGGCGATCCAGACTGAGGCGTGG-3'(SEQ.ID.NO.:28; antisense, containing the stop codon) and AP2 primer (Clontech) were used for the second round of PCR. DNA fragments generated by the 5' RACE-PCR were cloned into the pCRII-TOPO vector (Invitrogen) and sequenced using the SP6/T7 primers (Stratagene).
- Based on the sequence of the 5' RACE products, the full length RUP30 was cloned by RT-PCR, using primers:
 - 5'-GAGGTACAGCTGGCGATGCTGACAG-3' (SEQ.ID.NO.:29; sense, ATG as the initiation codon);
 - 5'-GTGGCCATGAGCCACCCTGAGCTCC-3' (SEQ.ID.NO.:30; antisense, 3' of the stop codon) and human pancreas Marathon-ReadyTM cDNA (Clontech) as template. Taq DNA polymerase (Stratagene) was used for the amplification in 50 μl reaction by the following cycle with step 2 to step 4 repeated 35 times: 94°C for 40 seconds; 94°C for 20 seconds; 64°C for 20 seconds; 72°C for 2 minutes; and 72°C for 5 minutes.
- A 1.2 Kb PCR fragment was isolated from a 1% agarose gel and cloned into the pCRII-TOPO vector (Invitrogen) and several clones were sequenced using the ABI Big Dye Terminator kit (P.E. Biosystems). See, SEQ.ID.NO.:5 for the nucleic acid sequence and SEQ.ID.NO.:6 for the putative amino acid sequence.

d. hRUP31 (Seq. Id. Nos. 7 & 8)

5

The disclosed human RUP31 was identified based upon the use of GenBank database information. While searching the database, a cDNA clone with Accession Number AL356214 was identified as a human genomic sequence from chromosome 10.

The full length RUP31 was cloned by RT-PCR using primers:

- 5'-GGAATGTCCACTGAATGCGCGCGG-3' (SEQ.ID.NO.:31; sense, containing the initiation codon),
- 5'-AGCTCGCCAGGTGTGAGAAACTCGG-3' (SEQ.ID.NO.:32; antisense, 3' of stop codon) and human mammary gland Marathon-ReadyTM cDNA (Clontech) as template.
- AdvantageTM cDNA polymerase (Clontech) was used for the amplification in 50 μl reaction by the following cycle with step 2 to step 4 repeated 35 times: 94°C for 40 sec; 94°C for 20 sec; 66°C for 20 sec; 72°C for 1 min 30 sec; and 72°C for 5 min.
 - A 1.1 kb PCR fragment was isolated from a 1% agarose gel and cloned into the pCRII-TOPO vector (Invitrogen) and sequenced using the ABI Big Dye Terminator Kit (P.E. Biosystems). See, SEQ.ID.NO.:7 for the nucleic acid sequence and SEQ.ID.NO.:8 for the putative amino acid sequence.

e. hRUP32 (Seq. Id. Nos. 9 & 10)

The disclosed human RUP32 was identified based upon the use of GenBank database information. While searching the database, a cDNA clone with Accession Number AL513524 was identified as a human genomic sequence from chromosome 6.

The full length RUP32 was cloned by PCR using primers:

- 5'-GCGTTATGAGCAGCAATTCATCCCTGCTGG-3' (SEQ.ID.NO.:33; sense, containing the initiation codon),
- 5'-GTATCCTGAACTTCGTCTATACAACTGC-3' (SEQ.ID.NO.:34; antisense)
- and human genomic DNA (Clontech) as template. TaqPlus[®] Precision DNA polymerase (Stratagene) was used for the amplification by the following cycle with step 2 to step 4 repeated 35 times: 94°C for 3 min; 94°C for 20 sec; 58°C for 20 sec; 72°C for 1 min 30 sec; and 72°C for 7 min.
 - A 1.06 kb PCR fragment was isolated from a 1% agarose gel and cloned into the pCRII-TOPO vector (Invitrogen) and sequenced using the ABI Big Dye Terminator Kit (P.E. Biosystems). See, SEQ.ID.NO.:9 for the nucleic acid sequence and SEQ.ID.NO.:10 for the putative amino acid sequence.

f. hRUP33 (Seq. Id. Nos. 11 & 12)

5

The disclosed human RUP33 was identified based upon the use of GenBank database information. While searching the database, a cDNA clone with Accession Number AL513524 was identified as a human genomic sequence from chromosome 6.

The full length RUP33 was cloned by PCR using primers:

- 5'-CCCTCAGGAATGATGCCCTTTTGCCACAA-3' (SEQ.ID.NO.:35; sense, containing the initiation codon),
- 5'-ATCCATGTGGTTGGTGCATGTGGTTCGT-3' (SEQ.ID.NO.:36; antisense) and human genomic DNA (Clontech) as template. TaqPlus® Precision DNA polymerase (Stratagene) was used for the amplification by the following cycle with step 2 to step 4 repeated 35 times: 94°C for 3 min; 94°C for 20 sec; 56°C for 20 sec; 72°C for 1 min 30 sec; and 72°C for 7 min.
- A 1.1 kb PCR fragment was isolated from a 1% agarose gel and cloned into the pCRII-TOPO vector (Invitrogen) and sequenced using the ABI Big Dye Terminator Kit (P.E. Biosystems). See, SEQ.ID.NO.:11 for the nucleic acid sequence and SEQ.ID.NO.:12 for the putative amino acid sequence.

g. hRUP34 (Seq. Id. Nos. 13 & 14)

The disclosed human RUP34 was identified based upon the use of GenBank database information. While searching the database, a cDNA clone with Accession Number AL513524 was identified as a human genomic sequence from chromosome 6.

The full length RUP34 was cloned by PCR using primers:

- 5'-AAACAACAACAGCAGCAGCATGACCAGC-3' (SEQ.ID.NO.:37; sense, containing the initiation codon),
- 5'-ACATAGAGACAAGTGACATGTGTGAACCAC-3' (SEQ.ID.NO.:38; antisense)
- and human genomic DNA (Clontech) as template. TaqPlus® Precision DNA polymerase (Stratagene) was used for the amplification by the following cycle with step 2 to step 4 repeated 35 times: 94°C for 3 min; 94°C for 20 sec; 60°C for 20 sec; 72°C for 1 min 30 sec; and 72°C for 7 min.
- A 1.27 kb PCR fragment was isolated from a 1% agarose gel and cloned into the pCRII-TOPO vector (Invitrogen) and sequenced using the ABI Big Dye Terminator Kit (P.E. Biosystems). See, SEQ.ID.NO.:13 for the nucleic acid sequence and SEQ.ID.NO.:14 for the putative amino acid sequence.

h. hRUP35 (Seq. Id. Nos. 15 & 16)

The disclosed human RUP35 was identified based upon the use of GenBank database information. While searching the database, a cDNA clone with Accession Number AC021089 was identified as a human genomic sequence from chromosome 16.

The 5' sequence of RUP35 was determined by 5' RACE- PCR with a human fetal brain Marathon-Ready™ cDNA (Clontech) as template. Oligonucleotide 5'-GGTATGAGACCGTGTGGTACTTGAGC-3' (SEQ.ID.NO.:39; sense) and AP1 primer (Clontech) were used for the first round of RT-PCR and oligonucleotide 5'-GTGGCAGACAGCGATATACCTGTCAATGG-3' (SEQ.ID.NO.:40; antisense)

and AP2 primer (Clontech) were used for the second round of PCR. DNA fragments generated by the 5' RACE-PCR were cloned into the pCRII-TOPO vector (Invitrogen) and sequenced using the SP6/T7 primers (Stratagene).

Based upon the sequence of the 5' RACE products, the full length RUP35 was cloned by RT-PCR, using primers

5'-GCGCTCATGGAGCACACGCACGCCCAC-3' (SEQ.ID.NO.:41; sense, ATG as the initiation codon) and

5'-GAGGCAGTAGTTGCCACACCTATGG-3' (SEQ.ID.NO.:42; antisense, 3' of the stop codon) and human brain Marathon-ReadyTM cDNA (Clontech) as template. AdvantageTM cDNA polymerase (Clontech) was used for the amplification in 100 μl reaction by the following cycle with step 2 to step 4 repeated 45 times: 95°C for 2 min; 95°C for 20 sec; 60°C for 20 sec; 72°C for 1 min 30 sec; and 72°C for 5 min.

A 1.0 kb PCR fragment was isolated from a 1% agarose gel and cloned into the pCRII-TOPO vector (Invitrogen) and sequenced using the ABI Big Dye Terminator Kit (P.E. Biosystems). See, SEQ.ID.NO.:15 for the nucleic acid sequence and SEQ.ID.NO.:16 for the putative amino acid sequence.

i. hRUP36 (Seq. Id. Nos. 17 & 18)

20

25

30

The disclosed human RUP36 was identified based upon the use of GenBank database information. While searching the database, a cDNA clone with Accession Number AC090099 was identified as a human genomic sequence from chromosome 11.

The full length RUP36 was cloned by PCR using primers:
5'- CATCTGGTTTGTGTTCCCAGGGGCACCAG -3' (SEQ.ID.NO.:43; sense, 5' of start codon).

5'-

15

25

GACAGTGTTGCTCTCAAAGTCCCGTCTGACTG -3' (SEQ.ID.NO.:44; antisense, 3' of stop codon) and human genomic DNA (Clontech) as template. TaqPlus[®] Precision DNA polymerase (Stratagene) was used for the amplification in a 50µl reaction by the following cycle with step 2 to step 4 repeated 30 times: 95°C, 5 min; 95°C for 30 sec; 70°C for 30 sec; 72°C for 1 min 30 sec; and 72°C for 7 min.

A 1.0 kb PCR fragment was isolated from a 1% agarose gel and cloned into the pCRII-TOPO vector (Invitrogen) and sequenced using the ABI Big Dye Terminator Kit (P.E. Biosystems). See, SEQ.ID.NO.:17 for the nucleic acid sequence and SEQ.ID.NO.:18 for the putative amino acid sequence.

j. hRUP37 (Seq. Id. Nos. 19 & 20)

The disclosed human RUP37 was identified based upon the use of GenBank database information. While searching the database, a cDNA clone with Accession Number AC090099 was identified as a human genomic sequence from chromosome 11.

The full length RUP37 was cloned by PCR using primers:

5'-CTGTTTCCAGGGTCATCAGACTGGG-3' (SEQ.ID.NO.:45; sense);

5'-GCAGCATTGCTCTCAAAGTCCTGTCTG-3' (SEQ.ID.NO.:46; antisense)

and human genomic DNA (Clontech) as template. TaqPlus[®] Precision DNA polymerase

(Stratagene) was used for the amplification by the following cycle with step 2 to step 4

repeated 35 times: 95°C for 5 min; 95°C for 30 sec; 62°C for 30 sec; 72°C for 1 min 30 sec;

and 72°C for 7 min.

A 969 base pair was isolated from a 1%% agarose gel and cloned into the pCRII-TOPO vector (Invitrogen) and sequenced using the ABI Big Dye Terminator Kit (P.E. Biosystems). See, SEQ.ID.NO.:19 for the nucleic acid sequence and SEQ.ID.NO.:20 for the putative amino acid sequence.

Example 2

PREPARATION OF NON-ENDOGENOUS, CONSTITUTIVELY ACTIVATED GPCRS

Those skilled in the art are credited with the ability to select techniques for mutation of a nucleic acid sequence. Presented below are approaches utilized to create non-endogenous versions of several of the human GPCRs disclosed above. The mutations disclosed below are based upon an algorithmic approach whereby the 16th amino acid (located in the IC3 region of the GPCR) from a conserved proline (or an

endogenous, conservative substitution therefore) residue (located in the TM6 region of the GPCR, near the TM6/IC3 interface) is mutated, preferably to an alanine, histimine, arginine or lysine amino acid residue, most preferably to a lysine amino acid residue.

1. Transformer Site-Directed ™ Mutagenesis

Preparation of non-endogenous human GPCRs may be accomplished on human GPCRs using, *inter alia*, Transformer Site-DirectedTM Mutagenesis Kit (Clontech) according to the manufacturer instructions. In some embodiments two mutagenesis primers are used, preferably a lysine mutagenesis oligonucleotide that creates the lysine mutation, and a selection marker oligonucleotide. For convenience, the codon mutation to be incorporated into the human GPCR is also noted, in standard form (Table D):

TABLE D

Receptor Identifier	Codon Mutation
hRUP28	V274K
hRUP29	T249K
hRUP30	R232K
hRUP31	M294K
hRUP32	F220K
hRUP34	A238K
hRUP35	Y215K
hRUP36	L294K
· hRUP37	T219K

15 Example 3

5

10

RECEPTOR EXPRESSION

Although a variety of cells are available to the art-skilled for the expression of proteins, it is preferred that mammalian cells be utilized. The primary reason for this is predicated upon practicalities, *i.e.*, utilization of, *e.g.*, yeast cells for the expression of a GPCR, while possible, introduces into the protocol a non-mammalian cell which may not (indeed, in the case of yeast, does not) include the receptor-coupling, genetic-mechanism

and secretary pathways that have evolved for mammalian systems — thus, results obtained in non-mammalian cells, while of potential use, are not as preferred as those obtained using mammalian cells. Of the mammalian cells, COS-7, 293 and 293T cells are particularly preferred, although the specific mammalian cell utilized can be predicated upon the particular needs of the artisan.

a. Transient Transfection

On day one, $6x10^6$ cells/10 cm dish of 293 cells well were plated out. On day two, two reaction tubes were prepared (the proportions to follow for each tube are per plate): tube A was prepared by mixing 4µg DNA (e.g., pCMV vector; pCMV vector with receptor cDNA, etc.) in 0.5 ml serum free DMEM (Gibco BRL); tube B was prepared by mixing 24µl lipofectamine (Gibco BRL) in 0.5ml serum free DMEM. Tubes A and B were admixed by inversion (several times), followed by incubation at room temperature for 30-45min. The admixture is referred to as the "transfection mixture". Plated 293 cells were washed with 1XPBS, followed by addition of 5 ml serum free DMEM. One ml of the transfection mixture were added to the cells, followed by incubation for 4hrs at 37°C/5% CO₂. The transfection mixture was removed by aspiration, followed by the addition of 10ml of DMEM/10% Fetal Bovine Serum. Cells were incubated at 37°C/5% CO₂. After 48hr incubation, cells were harvested and utilized for analysis.

b. Stable Cell Lines

20

30

Approximately 12x10⁶ 293 cells will be plated on a 15cm tissue culture plate, and grown in DME High Glucose Medium containing 10% fetal bovine serum and one percent sodium pyruvate, L-glutamine, and antibiotics. Twenty-four hours following plating of 293 cells (to approximately ~80% confluency), the cells will be transfected using 12μg of DNA. The 12μg of DNA is combined with 60μl of lipofectamine and 2mL of DME High Glucose Medium without serum. The medium will be aspirated from the plates and the cells washed once with medium without serum. The DNA, lipofectamine, and medium mixture will be added to the plate along with 10mL of medium without serum. Following incubation at 37°C for four to five hours, the medium will be aspirated and 25ml of medium containing serum will be added. Twenty-four hours following transfection, the medium will be aspirated again, and fresh medium with serum will be added. Forty-eight hours following transfection, the medium will be aspirated and medium with serum will be added containing geneticin (G418 drug) at a final concentration of 500μg/mL. The transfected cells will then

undergo selection for positively transfected cells containing the G418 resistant gene. The medium will be replaced every four to five days as selection occurs. During selection, cells will be grown to create stable pools, or split for stable clonal selection.

Example 4

10

15

20

25

30

5 ASSAYS FOR DETERMINATION OF CONSTITUTIVE ACTIVITY OF NON-ENDOGENOUS GPCRS

A variety of approaches are available for assessment of constitutive activity of the non-endogenous human GPCRs. The following are illustrative; those of ordinary skill in the art are credited with the ability to determine those techniques that are preferentially beneficial for the needs of the artisan.

1. Membrane Binding Assays: [35S]GTPγS Assay

When a G protein-coupled receptor is in its active state, either as a result of ligand binding or constitutive activation, the receptor couples to a G protein and stimulates the release of GDP and subsequent binding of GTP to the G protein. The alpha subunit of the G protein-receptor complex acts as a GTPase and slowly hydrolyzes the GTP to GDP, at which point the receptor normally is deactivated. Constitutively activated receptors continue to exchange GDP for GTP. The non-hydrolyzable GTP analog, [35S]GTPγS, can be utilized to demonstrate enhanced binding of [35S]GTPγS to membranes expressing constitutively activated receptors. Advantages of using [35S]GTPγS binding to measure constitutive activation include but are not limited to the following: (a) it is generically applicable to all G protein-coupled receptors; (b) it is proximal at the membrane surface making it less likely to pick-up molecules which affect the intracellular cascade.

The assay takes advantage of the ability of G protein coupled receptors to stimulate [³⁵S]GTPγS binding to membranes expressing the relevant receptors. The assay can, therefore, be used in the direct identification method to screen candidate compounds to constitutively activated G protein-coupled receptors. The assay is generic and has application to drug discovery at all G protein-coupled receptors.

The [35S]GTPγS assay is incubated in 20 mM HEPES and between 1 and about 20mM MgCl₂ (this amount can be adjusted for optimization of results, although 20mM is preferred) pH 7.4, binding buffer with between about 0.3 and about 1.2 nM [35S]GTPγS

(this amount can be adjusted for optimization of results, although 1.2 is preferred) and 12.5 to 75 μ g membrane protein (e.g., 293 cells expressing the G_s Fusion Protein; this amount can be adjusted for optimization) and 10 μ M GDP (this amount can be changed for optimization) for 1 hour. Wheatgerm agglutinin beads (25 μ l; Amersham) will then be added and the mixture incubated for another 30 minutes at room temperature. The tubes will be then centrifuged at 1500 x g for 5 minutes at room temperature and then counted in a scintillation counter.

2. Adenylyl Cyclase

10

15

20

25

A Flash PlateTM Adenylyl Cyclase kit (New England Nuclear; Cat. No. SMP004A) designed for cell-based assays can be modified for use with crude plasma membranes. The Flash Plate wells can contain a scintillant coating which also contains a specific antibody recognizing cAMP. The cAMP generated in the wells can be quantitated by a direct competition for binding of radioactive cAMP tracer to the cAMP antibody. The following serves as a brief protocol for the measurement of changes in cAMP levels in whole cells that express the receptors.

Transfected cells will be harvested approximately twenty four hours after transient transfection. Media will be carefully aspirated and discarded. Ten ml of PBS will gently be added to each dish of cells followed by careful aspiration. One ml of Sigma cell dissociation buffer and 3ml of PBS will be added to each plate. Cells will be pipetted off the plate and the cell suspension collected into a 50ml conical centrifuge tube. Cells will be centrifuged at room temperature at 1,100 rpm for 5 min. The cell pellet will be carefully resuspended into an appropriate volume of PBS (about 3ml/plate). The cells will be then counted using a hemocytometer and additional PBS will be added to give the appropriate number of cells (to a final volume of about 50µl/well).

cAMP standards and Detection Buffer (comprising 1 µCi of tracer [\$^{125}\$I cAMP (50 µl] to 11 ml Detection Buffer) will be prepared and maintained in accordance with the manufacturer's instructions. Assay Buffer will be prepared fresh for screening and contained 50µl of Stimulation Buffer, 3µl of test compound (12µM final assay concentration) and 50µl cells, Assay Buffer will be stored on ice until utilized. The assay will be initiated by addition of 50µl of cAMP standards to appropriate wells followed by addition of 50µl of PBSA to wells H-11 and H12. Fifty µl of Stimulation Buffer will be added to all wells. DMSO (or selected candidate compounds) will be added to appropriate

wells using a pin tool capable of dispensing 3μl of compound solution, with a final assay concentration of 12μM test compound and 100μl total assay volume. The cells will then be added to the wells and incubated for 60 min at room temperature. One hundred μl of Detection Mix containing tracer cAMP will then be added to the wells. Plates will be incubated for an additional 2 hours followed by counting in a Wallac MicroBetaTM scintillation counter. Values of cAMP/well will then be extrapolated from a standard cAMP curve which will be contained within each assay plate.

3. Cell-Based cAMP for G_i Coupled Target GPCRs

10

25

TSHR is a Gs coupled GPCR that causes the accumulation of cAMP upon activation. TSHR will be constitutively activated by mutating amino acid residue 623 (i.e., changing an alanine residue to an isoleucine residue). A Gi coupled receptor is expected to inhibit adenylyl cyclase, and, therefore, decrease the level of cAMP production, which can make assessment of cAMP levels challenging. An effective technique for measuring the decrease in production of cAMP as an indication of constitutive activation of a Gi coupled receptor can be accomplished by co-transfecting, most preferably, non-endogenous, constitutively activated TSHR (TSHR-A623I) (or an endogenous, constitutively active G_s coupled receptor) as a "signal enhancer" with a Gi linked target GPCR to establish a baseline level of cAMP. Upon creating a non-endogenous version of the Gi coupled receptor, this non-endogenous version of the target GPCR is then co-transfected with the signal enhancer, and it is this material that can be used for screening. This approach will be utilized to effectively generate a signal when a cAMP assay is used; this approach is preferably used in the direct identification of candidate compounds against Gi coupled receptors. It is noted that for a Gi coupled GPCR, when this approach is used, an inverse agonist of the target GPCR will increase the cAMP signal and an agonist will decrease the cAMP signal.

On day one, $2x10^4$ 293 cells/well will be plated out. On day two, two reaction tubes will be prepared (the proportions to follow for each tube are per plate): tube A will be prepared by mixing 2ug DNA of each receptor transfected into the mammalian cells, for a total of 4ug DNA (e.g., pCMV vector; pCMV vector with mutated THSR (TSHR-A623I); TSHR-A623I and GPCR, etc.) in 1.2ml serum free DMEM (Irvine Scientific, Irvine, CA); tube B will be prepared by mixing 120µl lipofectamine (Gibco BRL) in 1.2ml serum free

DMEM. Tubes A and B will then be admixed by inversion (several times), followed by incubation at room temperature for 30-45min. The admixture is referred to as the "transfection mixture". Plated 293 cells will be washed with 1XPBS, followed by addition of 10ml serum free DMEM. 2.4ml of the transfection mixture will then be added to the cells, followed by incubation for 4hrs at 37°C/5% CO₂. The transfection mixture will then be removed by aspiration, followed by the addition of 25ml of DMEM/10% Fetal Bovine Serum. Cells will then be incubated at 37°C/5% CO₂. After 24hr incubation, cells will be harvested and utilized for analysis.

A Flash Plate™ Adenylyl Cyclase kit (New England Nuclear; Cat. No. SMP004A) although designed for cell-based assays, can be modified for use with crude plasma membranes depending on the need of the skilled artisan. The Flash Plate wells will contain a scintillant coating which also contains a specific antibody recognizing cAMP. The cAMP generated in the wells can be quantified by a direct competition for binding of radioactive cAMP tracer to the cAMP antibody. The following serves as a brief protocol for the measurement of changes in cAMP levels in whole cells that express the receptors.

10

15

20

25

30

Transfected cells will be harvested approximately twenty four hours after transient transfection. Media will be carefully aspirated and discarded. Ten ml of PBS will be gently added to each dish of cells followed by careful aspiration. One ml of Sigma cell dissociation buffer and 3ml of PBS will be added to each plate. Cells will be pipetted off the plate and the cell suspension will be collected into a 50ml conical centrifuge tube. Cells will be centrifuged at room temperature at 1,100 rpm for 5 min. The cell pellet will be carefully re-suspended into an appropriate volume of PBS (about 3ml/plate). The cells will then be counted using a hemocytometer and additional PBS is added to give the appropriate number of cells (to a final volume of about 50µl/well).

cAMP standards and Detection Buffer (comprising 1 µCi of tracer [125] cAMP (50 µl] to 11 ml Detection Buffer) will be prepared and maintained in accordance with the manufacturer's instructions. Assay Buffer should be prepared fresh for screening and contained 50µl of Stimulation Buffer, 3µl of test compound (12µM final assay concentration) and 50µl cells, Assay Buffer can be stored on ice until utilized. The assay can be initiated by addition of 50µl of cAMP standards to appropriate wells followed by addition of 50µl of PBSA to wells H-11 and H12. Fifty µl of Stimulation Buffer will be added to all wells. Selected compounds (e.g., TSH) will be added to appropriate wells

using a pin tool capable of dispensing 3µl of compound solution, with a final assay concentration of 12µM test compound and 100µl total assay volume. The cells will then be added to the wells and incubated for 60 min at room temperature. One hundred µl of Detection Mix containing tracer cAMP will then be added to the wells. Plates will then be incubated additional 2 hours followed by counting in a Wallac MicroBeta scintillation counter. Values of cAMP/well will then be extrapolated from a standard cAMP curve which is contained within each assay plate.

4. Reporter-Based Assays

10

20

25

30

a. CRE-LUC Reporter Assay (G_s -associated receptors)

293 and 293T cells will be plated-out on 96 well plates at a density of 2 x 10⁴ cells per well and will be transfected using Lipofectamine Reagent (BRL) the following day according to manufacturer instructions. A DNA/lipid mixture will be prepared for each 6-well transfection as follows: 260ng of plasmid DNA in 100µl of DMEM are gently mixed with 2µl of lipid in 100µl of DMEM (the 260ng of plasmid DNA consisted of 200ng of a 8xCRE-Luc reporter plasmid, 50ng of pCMV comprising endogenous receptor or non-endogenous receptor or pCMV alone, and 10ng of a GPRS expression plasmid (GPRS in pcDNA3 (Invitrogen)). The 8XCRE-Luc reporter plasmid is prepared as follows: vector SRIF-β-gal will be obtained by cloning the rat somatostatin promoter (-71/+51) at BglV-HindIII site in the pβgal-Basic Vector (Clontech). Eight (8) copies of cAMP response element will be obtained by PCR from an adenovirus template AdpCF126CCRE8 (see, 7 Human Gene Therapy 1883 (1996)) and cloned into the SRIFβ-gal vector at the Kpn-BglV site, resulting in the 8xCRE-β-gal reporter vector. The 8xCRE-Luc reporter plasmid will be generated by replacing the beta-galactosidase gene in the 8xCRE-β-gal reporter vector with the luciferase gene obtained from the pGL3basic vector (Promega) at the HindIII-BamHI site. Following 30 min. incubation at room temperature, the DNA/lipid mixture will be diluted with 400 µl of DMEM and 100µl of the diluted mixture will be added to each well. One hundred µl of DMEM with 10% FCS will be added to each well after a 4hr incubation in a cell culture incubator. The following day the transfected cells will be changed with 200 µl/well of DMEM with 10% FCS. Eight hours later, the wells will be changed to 100 µl /well of DMEM without phenol red, after one wash with PBS. Luciferase activity will be measured the next day using the LucLite™ reporter gene assay kit (Packard) following manufacturer's

instructions and read on a 1450 MicroBeta™ scintillation and luminescence counter (Wallac).

b. AP1 reporter assay (G_q-associated receptors)

A method to detect G_q stimulation depends on the known property of G_{q^-} dependent phospholipase C to cause the activation of genes containing AP1 elements in their promoter. A PathdetectTM AP-1 cis-Reporting System (Stratagene, Catalogue # 219073) can be utilized following the protocol set forth above with respect to the CREB reporter assay, except that the components of the calcium phosphate precipitate were 410 ng pAP1-Luc, 80 ng pCMV-receptor expression plasmid, and 20 ng CMV-SEAP.

c. SRF-Luc Reporter Assay (G_q- associated receptors)

One method to detect G_q stimulation depends on the known property of G_qdependent phospholipase C to cause the activation of genes containing serum response factors in their promoter. A Pathdetect™ SRF-Luc-Reporting System (Stratagene) can be utilized to assay for Gq coupled activity in, e.g., COS7 cells. Cells are transfected with the plasmid components of the system and the indicated expression plasmid encoding endogenous or non-endogenous GPCR using a Mammalian Transfection™ Kit (Stratagene, Catalogue #200285) according to the manufacturer's instructions. Briefly, 410 ng SRF-Luc, 80 ng pCMV-receptor expression plasmid and 20 ng CMV-SEAP (secreted alkaline phosphatase expression plasmid; alkaline phosphatase activity is measured in the media of transfected cells to control for variations in transfection efficiency between samples) are combined in a calcium phosphate precipitate as per the manufacturer's instructions. Half of the precipitate is equally distributed between 3 wells in a 96-well plate, kept on the cells in a serum free media for 24 hours. The last 5 hours the cells are incubated with 1µM Angiotensin, where indicated. Cells are then lysed and assayed for luciferase activity using a Luclite™ Kit (Packard, Cat. # 6016911) and "Trilux 1450 Microbeta" liquid scintillation and luminescence counter (Wallac) as per the manufacturer's instructions. The data can be analyzed using GraphPad Prism™ 2.0a (GraphPad Software Inc.).

15

20

25

30

d. Intracellular IP₃ Accumulation Assay (G₀-associated receptors)

On day 1, cells comprising the receptors (endogenous and/or nonendogenous) are plated onto 24 well plates, usually 1x10⁵ cells/well (although his number can be optimized. On day 2 cells are transfected by firstly mixing 0.25ug DNA in 50 µl serum free DMEM/well and 2 µl lipofectamine in 50 µl serum free DMEM/well. The solutions are gently mixed and incubated for 15-30 min at room temperature. Cells are then washed with 0.5 ml PBS and 400 µl of serum free media and then mixed with the transfection media and added to the cells. The cells are incubated for 3-4 hrs at 37°C/5%CO2 and then the transfection media is removed and replaced with 1ml/well of regular growth media. On day 3 the cells are labeled with ³H-myo-inositol. Briefly, the media is removed and the cells are washed with 0.5 ml PBS. Then 0.5 ml inositolfree/serum free media (GIBCO BRL) are added/well with 0.25 µCi of ³H-myo-inositol/ well and the cells incubated for 16-18 hrs overnight at 37°C/5%CO₂. On Day 4 the cells are washed with 0.5 ml PBS and 0.45 ml of assay medium is added containing inositolfree/serum free media 10 µM pargyline 10 mM lithium chloride or 0.4 ml of assay medium and 50 µl of 10x ketanserin (ket) to final concentration of 10µM. The cells are then incubated for 30 min at 37°C. The cells are then washed with 0.5 ml PBS and 200 µl of fresh/ice cold stop solution (1M KOH; 18 mM Na-borate; 3.8 mM EDTA) is added to each well. The solution is kept on ice for 5-10 min (or until cells are lysed) and then neutralized by 200 µl of fresh/ice cold neutralization solution (7.5 % HCL). The lysate is then transferred into 1.5 ml Eppendorf tubes and 1 ml of chloroform/methanol (1:2) is added/tube. The solution is vortexed for 15 sec and the upper phase is applied to a Biorad AG1-X8TM anion exchange resin (100-200 mesh). First, the resin is washed with water at 1:1.25 W/V and 0.9 ml of upper phase is loaded onto the column. The column is then washed with 10 ml of 5 mM myo-inositol and 10 ml of 5 mM Na-borate/60mM Naformate. The inositol tris phosphates are eluted into scintillation vials containing 10 ml of scintillation cocktail with 2 ml of 0.1 M formic acid/ 1 M ammonium formate. The columns are regenerated by washing with 10 ml of 0.1 M formic acid/3M ammonium formate and rinsed twice with dd H₂O and stored at 4°C in water.

20

25

. 30

Reference is made to Figure 1. In Figure 1, 293 cells were transfected with G_q protein containing a six amino acid deletion, " $G_q(del)$ "; G_q protein fused to a G_i protein, " $G_q(del)/G_i$ "; endogenous RUP32; and RUP32 with $G_q(del)$ ("RUP32+ $G_q(del)/G_i$ "). The data indicate, based upon measuring IP₃ accumulation of RUP32 co-transfection of

 $G_q(del)/G_i$, that RUP32 does not endogenously couple to G_q protein. However when RUP32 was co-transfected with $G_q(del)/G_i$ fusion protein, RUP32 was forced to couple to G_q protein. RUP27+ $G_q(del)/G_i$ evidence about a nine (9) fold increase in IP3 accumulation when compared to endogenous RUP32. This data demonstrates that the $G_q(del)/G_i$ Fusion Construct can be co-transfected with a GPCR and used to screen for agonists or inverse agonists.

Reference is made to Figure 2. In Figure 2, 293 cells were transfected with RUP35 and RUP36 receptor and compared to the control, pCMV. The data indicate that both RUP35 and RUP36 receptor are endogenously, constitutively active. RUP35 evidences about a six (6) fold increase in intracellular inositol phosphate accumulation when compared to pCMV and RUP36 evidences about a four (4) fold increase when compared to pCMV.

Example 5

5

15

20

25

30

FUSION PROTEIN PREPARATION

a. GPCR: G_s Fusion Construct

The design of the constitutively activated GPCR-G protein fusion construct can be accomplished as follows: both the 5' and 3' ends of the rat G protein G_sα (long form; Itoh, H. et al., 83 PNAS 3776 (1986)) is engineered to include a HindIII (5'-AAGCTT-3') sequence thereon. Following confirmation of the correct sequence (including the flanking HindIII sequences), the entire sequence is shuttled into pcDNA3.1(-) (Invitrogen, cat. no. V795-20) by subcloning using the HindIII restriction site of that vector. The correct orientation for the G_s\alpha sequence will be determined after subcloning into pcDNA3.1(-). The modified pcDNA3.1(-) containing the rat G_sα gene at HindIII sequence is then verified; this vector will then be available as a "universal" G_s a protein vector. The pcDNA3.1(-) vector contains a variety of well-known restriction sites upstream of the HindIII site, thus beneficially providing the ability to insert, upstream of the Gs protein, the coding sequence of an endogenous, constitutively active GPCR. This same approach can be utilized to create other "universal" G protein vectors, and, of course, other commercially available or proprietary vectors known to the artisan can be utilized. In some embodiments, the important criteria is that the sequence for the GPCR be upstream and in-frame with that of the G protein.

5

15

30

Spacers in the restriction sites between the G protein and the GPCR are optional. The sense and anti-sense primers included the restriction sites for XbaI and EcoRV, respectively, such that spacers (attributed to the restriction sites) exist between the G protein and the GPCR.

PCR will then be utilized to secure the respective receptor sequences for fusion within the G_sα universal vector disclosed above, using the following protocol for each: 100ng cDNA for GPCR will be added to separate tubes containing 2μl of each primer (sense and anti-sense), 3μl of 10mM dNTPs, 10μl of 10XTaqPlusTM Precision buffer, 1μl of TaqPlusTM Precision polymerase (Stratagene: #600211), and 80μl of water. Reaction temperatures and cycle times for the GPCR will be as follows with cycle steps 2 through 4 were repeated 35 times: 94°C for 1 min; 94°C for 30 seconds; 62°C for 20 sec; 72°C 1 min 40sec; and 72°C 5 min. PCR products will be run on a 1% agarose gel and then purified. The purified products will be digested with XbaI and EcoRV and the desired inserts purified and ligated into the G_s universal vector at the respective restriction sites. The positive clones will be isolated following transformation and determined by restriction enzyme digestion; expression using 293 cells will be accomplished following the protocol set forth *infra*. Each positive clone for GPCR- G_s Fusion Protein will be sequenced to verify correctness.

b. G_q(6 amino acid deletion)/G_i Fusion Construct

The design of a G_q (del)/G_i fusion construct was accomplished as follows: the N-terminal six (6) amino acids (amino acids 2 through 7), having the sequence of TLESIM (SEQ.ID.NO.:47) Gαq-subunit was deleted and the C-terminal five (5) amino acids, having the sequence EYNLV (SEQ.ID.NO.:48) was replaced with the corresponding amino acids of the Gαi Protein, having the sequence DCGLF (SEQ.ID.NO.:49). This fusion construct was obtained by PCR using the following primers:

5'-gatcAAGCTTCCATGGCGTGCTGCCTGAGCGAGG-3' (SEQ.ID.NO.:50) and 5'-gatcGGATCCTTAGAACAGGCCGCAGTCCTTCAGGTTCAGCTGCAGGATGGTG-3' (SEQ.ID.NO.:51) and Plasmid 63313 which contains the mouse Gαq-wild type version with a hemagglutinin tag as template. Nucleotides in lower caps are included as spacers.

TaqPlus® Precision DNA polymerase (Stratagene) was utilized for the amplification by the following cycles, with steps 2 through 4 repeated 35 times: 95°C for 2 min; 95°C for 20 sec; 56°C for 20 sec; 72°C for 2 min; and 72°C for 7 min. The PCR

product will be cloned into a pCRII-TOPO vector (Invitrogen) and sequenced using the ABI Big Dye Terminator kit (P.E. Biosystems). Inserts from a TOPO clone containing the sequence of the fusion construct will be shuttled into the expression vector pcDNA3.1(+) at the HindIII/BamHI site by a 2 step cloning process.

5 Example 6

TISSUE DISTRIBUTION OF THE DISCLOSED HUMAN GPCRS: RT-PCR

RT-PCR was applied to confirm the expression and to determine the tissue distribution of several novel human GPCRs. Oligonucleotides utilized were GPCR-specific and the human multiple tissue cDNA panels (MTC, Clontech) as templates. Taq DNA polymerase (Stratagene) were utilized for the amplification in a 40μ l reaction according to the manufacturer's instructions. Twenty μ l of the reaction will be loaded on a 1.5% agarose gel to analyze the RT-PCR products. Table E, below, lists the receptors, the cycle conditions and the primers utilized, and also lists exemplary diseases/disorders linked to the receptors.

TABLE E

Receptor	Cycle	5' Primer	3' Primer	DNA	Tissue Expression
Identifier	Conditions	(SEQ.ID.NO.)	(SEQ.ID.NO.)	Fragment	
	Min ('), Sec (")				
	Cycles 2-4		:		
	repeated 35	·			
	times				
hRUP28	94°C for 5 min;	GTCCTCACT	CTGCGTCCAC	710bp	heart; kidney; liver; lung
	94°C for 30 sec;	GGTGGCCAT	CAGAGTCAC		and pancreas
	58°C for 30 sec,	GTACTCC (52)	GTCTCC (53)		
	72°C for 1 min,				
	and 72°C for 7				
	min				
hRUP29	94°C for 5 min;	CTTGGATGTT	GTTTGTGGCT	690bp	leukocyte and ovary
	94°C for 30 sec;	TGGGCTGCC	AACGGCACA		
	58°C for 30 sec,	CTTCTGC (54)	AAACACAAT		
	72°C for 1 min,		TCC (55)	_	

PATENT

AIREINO					PAIENI
	and 72°C for 7				
	min				
hRUP30	94°C for 2 min;	CTGCTCACG	GTGGCCATG	690bp	pancreas
	94°C for 15 sec;	GTTGACCGC	AGCCACCCT		
	58°C for 20 sec,	TACACTGC	GAGCTCC (57)		
	72°C for 1 min,	(56)			
	and 72°C for 10				
	min				
hRUP31	95°C for 4 min;	CTTCTTCTCC	CCAAATCA	516bp	colon, lung, pancreas,
	95°C for 1 min;	GACGTCAAG	GTGTGCAA		thymus; cerebral cortex,
	52°C for 30 sec,	G (58)	ATCG (59)		hippocampus of brain,
	72°C for 1 min,				and fat cells
	and 72°C for 7				
	min				
hRUP32	95°C for 4 min;	TGAATGGGT	CAACGGTCT	527bp	thymus
	95°C for 1 min;	CCTGTGTGA	GACAACCTC		
	52°C for 30 sec,	AA (60)	CT (61)		
	72°C for 1 min,				
	and 72°C for 7				
	min				
hRUP34	95°C for 4 min;	TTGCTGTGAT	CAGGAAGCC	534bp	peripheral blood
	95°C for 1 min;	GTGGCATTTT	CATAAAGGC		leukocyte ('PBL"),
	52°C for 30 sec,	G (62)	ATCAA (63)		prostate and kidney
	72°C for 1 min,				
	and 72°C for 7				
	min				
hRUP35	95°C for 4 min;	ACATCACCT	CCAGCATCTT	557bp	thalamus
	95°C for 1 min;	GCTTCCTGA	GATGCAGTG		
	52°C for 30 sec,	CC (64)	T (65)		
	72°C for 1 min,				
	and 72°C for 7				
	min				
hRUP37	95°C for 4 min;	CCATCTCCA	GCTGTTAAG	517bp	testis, cerebral cortex
	95°C for 1 min;	AAATCCTCA	AGCGGACAG		and hippocampus
	52°C for 30 sec,	GTC (66)	GAAA (67)		
	72°C for 1 min,				
L	l	L	L		L

AREN-0309	PATENT					
and 72°C for 7						
min						

Diseases and disorders related to receptors located in these tissues or regions include, but are not limited to, cardiac disorders and diseases (e.g. thrombosis, myocardial infarction; atherosclerosis; cardiomyopathies); kidney disease/disorders (e.g., renal failure; renal tubular acidosis; renal glycosuria; nephrogenic diabetes insipidus; cystinuria; polycystic kidney disease); eosinophilia; leukocytosis; leukopenia; ovarian cancer; sexual dysfunction; polycystic ovarian syndrome; pancreatitis and pancreatic cancer; irritable bowel syndrome; colon cancer; Crohn's disease; ulcerative colitis; diverticulitis; Chronic Obstructive Pulmonary Disease (COPD); Cystic Fibrosis; pneumonia; pulmonary hypertension; tuberculosis and lung cancer, Parkinson's disease; movement disorders and ataxias; learning and memory disorders; eating disorders (e.g., anorexia; bulimia, etc.); obesity; cancers; thymoma; myasthenia gravis; circulatory disorders; prostate cancer; prostatitis; kidney disease/disorders(e.g., renal failure; renal tubular acidosis; renal glycosuria; nephrogenic diabetes insipidus; cystinuria; polycystic kidney disease); sensorimotor processing and arousal disorders; obsessive-compulsive disorders; testicular cancer; priapism; prostatitis; hernia; endocrine disorders; sexual dysfunction; allergies; depression; psychotic disorders; migraine; reflux; schizophrenia; ulcers; bronchospasm; epilepsy; prostatic hypertrophy; anxiety; rhinitis; angina; and glaucoma. Accordingly, the methods of the present invention may also be useful in the diagnosis and/or treatment of these and other diseases and disorders.

Example 7

15

20

25

30

Protocol: Direct Identification of Inverse Agonists and Agonists

A. [³⁵S]GTPγS Assay

Although endogenous, constitutively active GPCRs have been used for the direct identification of candidate compounds as, e.g., inverse agonists, for reasons that are not altogether understood, intra-assay variation can become exacerbated. In some embodiments a GPCR Fusion Protein, as disclosed above, is also utilized with a non-endogenous, constitutively activated GPCR. When such a protein is used, intra-assay variation appears to be substantially stabilized, whereby an effective signal-to-noise ratio is

obtained. This has the beneficial result of allowing for a more robust identification of candidate compounds. Thus, in some embodiments it is preferred that for direct identification, a GPCR Fusion Protein be used and that when utilized, the following assay protocols be utilized.

1. Membrane Preparation

Membranes comprising the constitutively active orphan GPCR Fusion Protein of interest and for use in the direct identification of candidate compounds as inverse agonists or agonists are preferably prepared as follows:

a. Materials

5

10

15

20

25

30

"Membrane Scrape Buffer" is comprised of 20mM HEPES and 10mM EDTA, pH 7.4; "Membrane Wash Buffer" is comprised of 20 mM HEPES and 0.1 mM EDTA, pH 7.4; "Binding Buffer" is comprised of 20mM HEPES, 100 mM NaCl, and 10 mM MgCl₂, pH 7.4

b. Procedure

All materials will be kept on ice throughout the procedure. Firstly, the media will be aspirated from a confluent monolayer of cells, followed by rinse with 10ml cold PBS, followed by aspiration. Thereafter, 5ml of Membrane Scrape Buffer will be added to scrape cells; this will be followed by transfer of cellular extract into 50ml centrifuge tubes (centrifuged at 20,000 rpm for 17 minutes at 4°C). Thereafter, the supernatant will be aspirated and the pellet will be resuspended in 30ml Membrane Wash Buffer followed by centrifugation at 20,000 rpm for 17 minutes at 4°C. The supernatant will then be aspirated and the pellet resuspended in Binding Buffer. The resuspended pellet will then be homogenized using a Brinkman Polytron™ homogenizer (15-20 second bursts until the material is in suspension). This is referred to herein as "Membrane Protein".

2. Bradford Protein Assay

Following the homogenization, protein concentration of the membranes will be determined, for example, using the Bradford Protein Assay (protein can be diluted to about 1.5mg/ml, aliquoted and frozen (-80°C) for later use; when frozen, protocol for use will be as follows: on the day of the assay, frozen Membrane Protein is thawed at room temperature, followed by vortex and then homogenized with a Polytron at about 12 x 1,000 rpm for about 5-10 seconds; it was noted that for multiple preparations, the homogenizer is thoroughly cleaned between homogenization of different preparations).

a. Materials

Binding Buffer (as discussed above); Bradford Dye Reagent; Bradford Protein Standard will be utilized, following manufacturer instructions (Biorad, cat. no. 500-0006).

b. Procedure

5

15

20

25

30

Duplicate tubes will be prepared, one including the membrane, and one as a control "blank". Each contains 800µl Binding Buffer. Thereafter, 10µl of Bradford Protein Standard (1mg/ml) will be added to each tube, and 10µl of membrane Protein will then be added to just one tube (not the blank). Thereafter, 200µl of Bradford Dye Reagent will be added to each tube, followed by vortexing. After five minutes, the tubes will be re-vortexed and the material therein will be transferred to cuvettes. The cuvettes will then be read using a CECIL 3041 spectrophotometer, at wavelength 595.

3. Direct Identification Assay

a. Materials

GDP Buffer consisted of 37.5 ml Binding Buffer and 2mg GDP (Sigma, cat. no. G-7127), followed by a series of dilutions in Binding Buffer to obtain 0.2 µM GDP (final concentration of GDP in each well was 0.1 µM GDP); each well comprising a candidate compound, has a final volume of 200µl consisting of 100µl GDP Buffer (final concentration, 0.1µM GDP), 50µl Membrane Protein in Binding Buffer, and 50µl [35S]GTPγS (0.6 nM) in Binding Buffer (2.5 µl [35S]GTPγS per 10ml Binding Buffer).

b. Procedure

Candidate compounds will be preferably screened using a 96-well plate format (these can be frozen at -80°C). Membrane Protein (or membranes with expression vector excluding the GPCR Fusion Protein, as control), will be homogenized briefly until in suspension. Protein concentration will then be determined using, for example, the Bradford Protein Assay set forth above. Membrane Protein (and controls) will then be diluted to 0.25mg/ml in Binding Buffer (final assay concentration, 12.5μg/well). Thereafter, 100 μl GDP Buffer is added to each well of a Wallac ScintistripTM (Wallac). A 5μl pin-tool will then be used to transfer 5 μl of a candidate compound into such well (i.e., 5μl in total assay volume of 200 μl is a 1:40 ratio such that the final screening concentration of the candidate compound is 10μM). Again, to avoid contamination, after each transfer step the pin tool is rinsed in three reservoirs comprising water (1X), ethanol (1X) and water (2X) – excess

liquid is shaken from the tool after each rinse and the tool is dried with paper and Kim wipes. Thereafter, 50 μl of Membrane Protein will be added to each well (a control well comprising membranes without the GPCR Fusion Protein was also utilized), and preincubated for 5-10 minutes at room temperature. Thereafter, 50 μl of [35S]GTPγS (0.6 nM) in Binding Buffer will be added to each well, followed by incubation on a shaker for 60 minutes at room temperature (again, in this example, plates were covered with foil). The assay will be stopped by spinning the plates at 4000 RPM for 15 minutes at 22°C. The plates will then be aspirated with an 8 channel manifold and sealed with plate covers. The plates will then be read on a Wallac 1450 using setting "Prot. #37" (as per manufacturer's instructions).

B. Cyclic AMP Assay

10

20

25

30

Another assay approach to directly identify candidate compound will be accomplished utilizing a cyclase-based assay. In addition to direct identification, this assay approach can be utilized as an independent approach to provide confirmation of the results from the [35S]GTPγS approach as set forth above.

A modified Flash PlateTM Adenylyl Cyclase kit (New England Nuclear; Cat. No. SMP004A) will be preferably utilized for direct identification of candidate compounds as inverse agonists and agonists to GPCRs in accordance with the following protocol.

Transfected cells will be harvested approximately three days after transfection. Membranes will be prepared by homogenization of suspended cells in buffer containing 20mM HEPES, pH 7.4 and 10mM MgCl₂. Homogenization will be performed on ice using a Brinkman PolytronTM for approximately 10 seconds. The resulting homogenate will be centrifuged at 49,000 X g for 15 minutes at 4°C. The resulting pellet will then be resuspended in buffer containing 20mM HEPES, pH 7.4 and 0.1 mM EDTA, homogenized for 10 seconds, followed by centrifugation at 49,000 X g for 15 minutes at 4°C. The resulting pellet will then be stored at -80°C until utilized. On the day of direct identification screening, the membrane pellet will slowly be thawed at room temperature, resuspended in buffer containing 20mM HEPES, pH 7.4 and 10mM MgCl₂, to yield a final protein concentration of 0.60mg/ml (the resuspended membranes will be placed on ice until use).

cAMP standards and Detection Buffer (comprising 2 μ Ci of tracer [125 I cAMP (100 μ I] to 11 ml Detection Buffer) will be prepared and maintained in accordance with the

manufacturer's instructions. Assay Buffer will be prepared fresh for screening and contain 20mM HEPES, pH 7.4, 10mM MgCl₂, 20mM phosphocreatine (Sigma), 0.1 units/ml creatine phosphokinase (Sigma), 50 μM GTP (Sigma), and 0.2 mM ATP (Sigma); Assay Buffer will be stored on ice until utilized.

Candidate compounds identified as per above (if frozen, thawed at room temperature) will be added, preferably, to 96-well plate wells (3 μ l/well; 12 μ M final assay concentration), together with 40 μ l Membrane Protein (30 μ g/well) and 50 μ l of Assay Buffer. This admixture will be incubated for 30 minutes at room temperature, with gentle shaking.

Following the incubation, 100µl of Detection Buffer will be added to each well, followed by incubation for 2-24 hours. Plates will then be counted in a Wallac MicroBetaTM plate reader using "Prot. #31" (as per manufacturer instructions).

C. Melanophore Screening Assay

5

10

20

25

30

A method for identifying candidate agonists or inverse agonists for a GPCR can be preformed by introducing tests cells of a pigment cell line capable of dispersing or aggregating their pigment in response to a specific stimulus and expressing an exogenous clone coding for the GCPR. A stimulant, e.g., light, sets an initial state of pigment disposition wherein the pigment is aggregated within the test cells if activation of the GPCR induces pigment dispersion. However, stimulating the cell with a stimulant to set an initial state of pigment disposition wherein the pigment is dispersed if activation of the GPCR induces pigment aggregation. The tests cells are then contacted with chemical compounds, and it is determined whether the pigment disposition in the cells changed from the initial state of pigment disposition. Dispersion of pigments cells due to the candidate compound coupling to the GPCR will appear dark on a petri dish, while aggregation of pigments cells will appear light.

Materials and methods will be followed according to the disclosure of U.S. Patent Number 5,462,856 and U.S. Patent Number 6,051,386, each of which are incorporated by reference.

Although a variety of expression vectors are available to those in the art, for purposes of utilization for both the endogenous and non-endogenous human GPCRs, in some embodiments it is preferred that the vector utilized be pCMV. This vector was deposited with the American Type Culture Collection (ATCC) on October 13, 1998 (10801)

University Blvd., Manassas, VA 20110-2209 USA) under the provisions of the Budapest Treaty for the International Recognition of the Deposit of Microorganisms for the Purpose of Patent Procedure. The DNA was tested by the ATCC and determined to be viable. The ATCC has assigned the following deposit number to pCMV: ATCC #203351.

References cited throughout this patent document, including co-pending and related patent applications, unless otherwise indicated, are fully incorporated herein by reference. Modifications and extension of the disclosed inventions that are within the purview of the skilled artisan are encompassed within the above disclosure and the claims that follow.

5

CLAIMS

What is claimed is:

5

25

 A G protein-coupled receptor encoded by an amino acid sequence of SEQ.ID.NO.:2.

- 2. A non-endogenous, constitutively activated version of the G protein-coupled receptor of claim 1.
- 3. A plasmid comprising a vector and the cDNA of SEQ.ID.NO.:1.
- 4. A host cell comprising the plasmid of claim 3.
- A G protein-coupled receptor encoded by an amino acid sequence of SEQ.ID.NO.:4.
 - 6. A non-endogenous, constitutively activated version of the G protein-coupled receptor of claim 5.
 - 7. A plasmid comprising a vector and the cDNA of SEQ.ID.NO.:3.
- 8. A host cell comprising the plasmid of claim 7.
 - A G protein-coupled receptor encoded by an amino acid sequence of SEQ.ID.NO.:6.
 - 10. A non-endogenous, constitutively activated version of the G protein-coupled receptor of claim 9.
- 20 11. A plasmid comprising a vector and the cDNA of SEQ.ID.NO.:5.
 - 12. A host cell comprising the plasmid of claim 11.
 - 13. A G protein-coupled receptor encoded by an amino acid sequence of SEQ.ID.NO.:8.
 - 14. A non-endogenous, constitutively activated version of the G protein-coupled receptor of claim 13.
 - 15. A plasmid comprising a vector and the cDNA of SEQ.ID.NO.:7.
 - 16. A host cell comprising the plasmid of claim 15.
 - 17. A G protein-coupled receptor encoded by an amino acid sequence of SEQ.ID.NO.:10.
- 30 18. A non-endogenous, constitutively activated version of the G protein-coupled receptor of claim 17.
 - 19. A plasmid comprising a vector and the cDNA of SEQ.ID.NO.:9.

- 20. A host cell comprising the plasmid of claim 19.
- 21. A G protein-coupled receptor encoded by an amino acid sequence of SEQ.ID.NO.:12.
- 22. A non-endogenous, constitutively activated version of the G protein-coupled receptor of claim 21.
- 23. A plasmid comprising a vector and the cDNA of SEQ.ID.NO.:11.
- 24. A host cell comprising the plasmid of claim 23.

5

15

- 25. A G protein-coupled receptor encoded by an amino acid sequence of SEQ.ID.NO.:14.
- 10 26. A non-endogenous, constitutively activated version of the G protein-coupled receptor of claim 25.
 - 27. A plasmid comprising a vector and the cDNA of SEQ.ID.NO.:13.
 - 28. A host cell comprising the plasmid of claim 27.
 - 29. A G protein-coupled receptor encoded by an amino acid sequence of SEQ.ID.NO.:16.
 - 30. A non-endogenous, constitutively activated version of the G protein-coupled receptor of claim 29.
 - 31. A plasmid comprising a vector and the cDNA of SEQ.ID.NO.:15.
 - 32. A host cell comprising the plasmid of claim 31.
- 33. A G protein-coupled receptor encoded by an amino acid sequence of SEQ.ID.NO.:18.
 - 34. A non-endogenous, constitutively activated version of the G protein-coupled receptor of claim 33.
 - 35. A plasmid comprising a vector and the cDNA of SEQ.ID.NO.:17.
- 25 36. A host cell comprising the plasmid of claim 35.
 - 37. A G protein-coupled receptor encoded by an amino acid sequence of SEQ.ID.NO.:20.
 - 38. A non-endogenous, constitutively activated version of the G protein-coupled receptor of claim 37.
- 39. A plasmid comprising a vector and the cDNA of SE.ID.NO.:19.
 - 40. A host cell comprising the plasmid of claim 39.

Figure 1



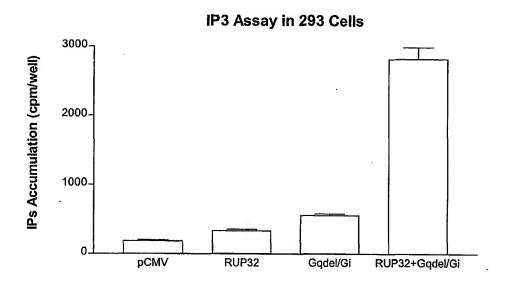


Figure 2

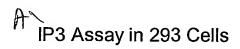

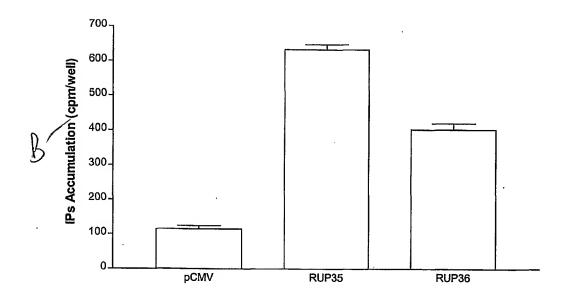

SEQUENCE LISTING

•

<110>	Arena Pharmaceuticals, Inc.
<120>	Endogenous And Non-Endogenous, Constitutively Activated Human G Protein-Coupled Receptors
<130>	AREN-0309
<150>	09/170,496
<151>	1998-10-13
<150>	PCT/US99/23938
<151>	1998-10-13
<150>	60/253,404
<151>	2000-11-27
<150>	60/255,366
<151>	2000-12-12
<150>	60/270,286
<151>	2001-02-20
<150>	60/282,365
<151>	2001-04-06
<150>	60/270,266
<151>	2001-02-20
<150>	60/282,032
<151>	2001-04-06
<150>	60/282,358
<151>	2001-04-06
<150>	60/282,356
<151>	2001-04-06
<150>.	60/290,917
<151>	2001-05-14
<150>	60/309,208
<151>	2001-07-31
<160>	67
<170>	PatentIn version 3.1
<212>	1 1002 DNA Artificial Sequence
<220> <223>	Novel Sequence
<400>	1
atgtgga	gct gcagctggtt caacggcaca gggctggtgg aggagctgcc tgcctgccag 60
gacctgc	agc tggggctgtc actgttgtcg ctgctgggcc tggtggtggg cgtgccagtg 120
ggcctgt	gct acaacgccct gctggtgctg gccaacctac acagcaaggc cagcatgacc 180
atgccgg	acg tgtactttgt caacatggca gtggcaggcc tggtgctcag cgccctggcc 240

cctgtgcacc tgctcggccc	cccgagctcc	cggtgggcgc	tgtggagtgt	gggcggcgaa	300
gtccacgtgg cactgcagat	ccccttcaat	gtgtcctcac	tggtggccat	gtactccacc	360
gccctgctga gcctcgacca	ctacatcgag	cgtgcactgc	cgcggaccta	catggccagc	420
gtgtacaaca cgcggcacgt	gtgcggcttc	gtgtggggtg	gcgcgctgct	gaccagcttc	480
tectegetge tettetacat	ctgcagccat	gtgtccaccc	gcgcgctaga	gtgcgccaag	540
atgcagaacg cagaagctgc	cgacgccacg	ctggtgttca	tcggctacgt	ggtgccagca	600
ctggccaccc tctacgcgct	ggtgctactc	tecegegtee	gcagggagga	cacgcccctg	660
gaccgggaca cgggccggct	ggagccctcg	gcacacaggc	tgctggtggc	caccgtgtgc	720
acgcagtttg ggctctggac	gccacactat	ctgatcctgc	tggggcacac	ggtcatcatc	780
tcgcgaggga agcccgtgga	cgcacactac	ctggggctac	tgcactttgt	gaaggatttc	840
tccaaactcc tggccttctc	cagcagcttt	gtgacaccac	ttctctaccg	ctacatgaac	900
cagagettee ecageaaget	ccaacggctg	atgaaaaagc	tgccctgcgg	ggaccggcac	960
tgctccccgg accacatggg	ggtgcagcag	gtgctggcgt	ag .		1002

<210> 2

<211> 333 <212> PRT

<213> Artificial Sequence

<220>

<223> Novel Sequence

<400> 2

Met Trp Ser Cys Ser Trp Phe Asn Gly Thr Gly Leu Val Glu Glu Leu

Pro Ala Cys Gln Asp Leu Gln Leu Gly Leu Ser Leu Leu Ser Leu Leu 20 25 30

Gly Leu Val Val Gly Val Pro Val Gly Leu Cys Tyr Asn Ala Leu Leu 35 40 45

Val Leu Ala Asn Leu His Ser Lys Ala Ser Met Thr Met Pro Asp Val

Tyr Phe Val Asn Met Ala Val Ala Gly Leu Val Leu Ser Ala Leu Ala 65 70 75 80

Pro Val His Leu Leu Gly Pro Pro Ser Ser Arg Trp Ala Leu Trp Ser

Val Gly Gly Glu Val His Val Ala Leu Gln Ile Pro Phe Asn Val Ser

Ser Leu Val Ala Met Tyr Ser Thr Ala Leu Leu Ser Leu Asp His Tyr Page 2

115 120 125

Ile	Glu 130	Arg	Ala	Leu	Pro	Arg 135	Thr	Tyr	Met	Ala	Ser 140	Val	Tyr	Asn	Thr	
Arg 145	His	Val	Cys	Gly	Phe 150	Val	Trp	Gly	Gly	Ala 155	Leu	Leu	Thr	Ser	Phe 160	
Ser	Ser	Leu	Leu	Phe 165	Tyr	Ile	Cys	Ser	His 170	Val	Ser	Thr	Arg	Ala 175	Leu	
Glu	Суз	Ala	Lys 180	Met	Gln	Asn	Ala	Glu 185	Ala	Ala	Asp	Ala	Thr 190	Leu	Val	
Phe	Ile	Gly 195	Tyr	Val	Val	Pro	Ala 200	Leu	Ala	Thr	Leu	Tyr 205	Ala	Leu	Val	
Leu	Leu 210	Ser	Arg	Val	Arg	Arg 215	Glu	Asp	Ťhr	Pro	Leu 220	Asp	Arg	Asp	Thr	
Gly 225	Arg	Leu	Glu	Pro	Ser 230	Ala	His	Arg	Leu	Leu 235	Val	Ala	Thr	Val	Cys 240	
Thr	Gln	Phe	Gly	Leu 245	Trp	Thr	Pro	His	Tyr 250	Leu	Ile	Leu	Leu	Gly 255	His	
Thr	Val	Ile	Ile 260	Ser	Arg	Gly	Lys	Pro 265	Val	Asp	Ala	His	Tyr 270	Leu	Gly	
Leu	Leu	His 275	Phe	Val	Lys	Asp	Phe 280	Ser	Lys	Leu	Leu	Ala 285	Phe	Ser	Ser	
Ser	Phe 290	Val	Thr	Pro	Leu	Leu 295	Tyr	Arg	Tyr	Met	Asn 300	Gln	Ser	Phe	Pro	
Ser 305	Lys	Leu	Gln	Arg	Leu 310	Met	Lys	Lys	Leu	Pro 315	Cys	Gly	Asp	Arg	His 320	
Cys	Ser	Pro	Asp	His 325	Met	Gly	Val	Gln	Gln 330	Val	Leu	Ala				
<210 <211 <212 <213	.> 9 !> 0	3 18 NA Iomo	sapi	.ens												
<400 atgc			caat	acct	c ca	ıggaa	ttcc	tct:	tgcg	atc	ctat	agtg	ac a	cccc	actta	60
atca	gcct	ct a	cttc	atag	rt go	ttat	tggc	ggg	ctgg	tgg	gtgt	catt	tc c	atto	ttttc	120
ctcc	atcagcctct acttcatagt gcttattggc gggctggtgg gtgtcatttc cattctttc 120 ctcctggtga aaatgaacac ccggtcagtg accaccatgg cggtcattaa cttggtggtg 180 Page 3															

gtccacagcg tttttctgct gacagtgcca tttcgcttga cctacctcat caagaagact 240 tggatgtttg ggctgccctt ctgcaaattt gtgagtgcca tgctgcacat ccacatgtac 300 ctcacgttcc tattctatgt ggtgatcctg gtcaccagat acctcatctt cttcaagtgc 360 aaagacaaag tggaattcta cagaaaactg catgctgtgg ctgccagtgc tggcatgtgg 420 acgctggtga ttgtcattgt ggtacccctg gttgtctccc ggtatggaat ccatgaggaa 480 tacaatgagg agcactgttt taaatttcac aaagagcttg cttacacata tgtgaaaatc 540 atcaactata tgatagtcat ttttgtcata gccgttgctg tgattctgtt ggtcttccag 600 gtcttcatca ttatgttgat ggtgcagaag ctacgccact ctttactatc ccaccaggag 660 ttctgggctc agctgaaaaa cctattttt ataggggtca tccttgtttg tttccttccc 720 taccagttct ttaggatcta ttacttgaat gttgtgacgc attccaatgc ctgtaacagc 780 aaggttgcat tttataacga aatcttcttg agtgtaacag caattagctg ctatgatttg 840 cttctctttg tctttggggg aagccattgg tttaagcaaa agataattgg cttatggaat 900 tgtgttttgt gccgttag 918

<210> 4

<211> 305

<212> PRT

<213> Homo sapiens

<400> 4

Met Pro Gly His Asn Thr Ser Arg Asn Ser Ser Cys Asp Pro Ile Val 1 5 15

Thr Pro His Leu Ile Ser Leu Tyr Phe Ile Val Leu Ile Gly Gly Leu 20 25 30

Val Gly Val Ile Ser Ile Leu Phe Leu Leu Val Lys Met Asn Thr Arg 35 40 45

Ser Val Thr Thr Met Ala Val Ile Asn Leu Val Val Val His Ser Val 50 60

Phe Leu Leu Thr Val Pro Phe Arg Leu Thr Tyr Leu Ile Lys Lys Thr 65 70 75 80

Trp Met Phe Gly Leu Pro Phe Cys Lys Phe Val Ser Ala Met Leu His 85 90 95

Ile His Met Tyr Leu Thr Phe Leu Phe Tyr Val Val Ile Leu Val Thr 100 105 110

Arg Tyr Leu Ile Phe Phe Lys Cys Lys Asp Lys Val Glu Phe Tyr Arg 115 120 125

Lys Leu His Ala Val Ala Ala Ser Ala Gly Met Trp Thr Leu Val Ile 130 135 140									
Val Ile Val Val Pro Leu Val Val Ser Arg Tyr Gly Ile His Glu Glu 145 150 155 160									
Tyr Asn Glu Glu His Cys Phe Lys Phe His Lys Glu Leu Ala Tyr Thr 165 170 175									
Tyr Val Lys Ile Ile Asn Tyr Met Ile Val Ile Phe Val Ile Ala Val 180 185 190									
Ala Val Ile Leu Leu Val Phe Gln Val Phe Ile Ile Met Leu Met Val 195 200 205									
Gln Lys Leu Arg His Ser Leu Leu Ser His Gln Glu Phe Trp Ala Gln 210 215 220									
Leu Lys Asn Leu Phe Phe Ile Gly Val Ile Leu Val Cys Phe Leu Pro 225 230 235 240									
Tyr Gln Phe Phe Arg Ile Tyr Tyr Leu Asn Val Val Thr His Ser Asn 245 250 255									
Ala Cys Asn Ser Lys Val Ala Phe Tyr Asn Glu Ile Phe Leu Ser Val 260 265 270									
Thr Ala Ile Ser Cys Tyr Asp Leu Leu Leu Phe Val Phe Gly Gly Ser 275 280 285									
His Trp Phe Lys Gln Lys Ile Ile Gly Leu Trp Asn Cys Val Leu Cys 290 295 300									
Arg 305 .									
<210> 5 <211> 1125 <212> DNA <213> Artificial Sequence									
<220> <223> Novel Sequence									
<400> 5 atgctgacag ggagctgcgg ggaccctcag aaaaagccac aggtgaccca ggactcaggg 60									
ccccagagca tggggcttga gggacgagag acagctggcc agccacgagt gaccctgctg 120									
cccacgccca acgtcagcgg gctgagccag gagtttgaaa gccactggcc agagatcgca 180									
gagaggtccc cgtgtgtggc tggcgtcatc cctgtcatct actacagtgt cctgctgggc 240									
Dage 5									

ttggggctgc ctgtcagcct cctgaccgca gtggccctgg cgcgccttgc caccaggacc 300 aggaggeeet ectaetaeta eettetggeg eteacageet eggatateat eateeaggtg 360 gtcatcgtgt tcgcgggctt cctcctgcag ggagcagtgc tggcccgcca ggtqccccaq 420 gctgtggtgc gcacggccaa catcctggag tttgctgcca accacgcctc agtctggatc 480 gccatcctgc tcacggttga ccgctacact gccctgtgcc acccctgca ccatcgggcc 540 gcctcgtccc caggccggac ccgccgggcc attgctgctg tcctgagtgc tgccctgttg 600 accggcatcc ccttctactg gtggctggac atgtggagag acaccgactc acccagaaca 660 etggacgagg tecteaagtg ggeteactgt etcaetgtet attteatece ttgtggegtg 720 ttcctggtca ccaactcggc catcatccac cggctacgga ggaggggccg gagtgggctg 780 cagccccggg tgggcaagag cacagccatc ctcctgggca tcaccacact gttcaccctc 840 ctgtgggcgc cccgggtctt cgtcatgctc taccacatgt acgtggcccc tgtccaccqq 900 gactggaggg tccacctggc cttggatgtg gccaatatgg tggccatgct ccacacggca 960 gccaacttcg gcctctactg ctttgtcagc aagactttcc gggccactgt ccgacaggtc 1020 atccacgatg cctacctgcc ctgcactttg gcatcacagc cagagggcat ggcggcgaag 1080 cctgtgatgg agcctccggg actccccaca ggggcagaag tgtag 1125

```
<210> 6
```

Met Leu Thr Gly Ser Cys Gly Asp Pro Gln Lys Lys Pro Gln Val Thr 1 5 10 15

Gln Asp Ser Gly Pro Gln Ser Met Gly Leu Glu Gly Arg Glu Thr Ala 20 25 30

Gly Gln Pro Arg Val Thr Leu Leu Pro Thr Pro Asn Val Ser Gly Leu 35 40 45

Ser Gln Glu Phe Glu Ser His Trp Pro Glu Ile Ala Glu Arg Ser Pro 50 55

Cys Val Ala Gly Val Ile Pro Val Ile Tyr Tyr Ser Val Leu Leu Gly 65 70 75 80

Leu Gly Leu Pro Val Ser Leu Leu Thr Ala Val Ala Leu Ala Arg Leu 85 90 95

Ala Thr Arg Thr Arg Pro Ser Tyr Tyr Tyr Leu Leu Ala Leu Thr
Page 6

<211> 374

<212> PRT

<213> Artificial Sequence

<220>

<223> Novel Sequence

<400> 6

Ala Ser Asp Ile Ile Ile Gln Val Val Ile Val Phe Ala Gly Phe Leu 115 Leu Gln Gly Ala Val Leu Ala Arg Gln Val Pro Gln Ala Val Val Arg Thr Ala Asn Ile Leu Glu Phe Ala Ala Asn His Ala Ser Val Trp Ile Ala Ile Leu Leu Thr Val Asp Arg Tyr Thr Ala Leu Cys His Pro Leu His His Arg Ala Ala Ser Ser Pro Gly Arg Thr Arg Arg Ala Ile Ala 180 185 190 Ala Val Leu Ser Ala Ala Leu Leu Thr Gly Ile Pro Phe Tyr Trp Trp Leu Asp Met Trp Arg Asp Thr Asp Ser Pro Arg Thr Leu Asp Glu Val Leu Lys Trp Ala His Cys Leu Thr Val Tyr Phe Ile Pro Cys Gly Val Phe Leu Val Thr Asn Ser Ala Ile Ile His Arg Leu Arg Arg Gly Arg Ser Gly Leu Gln Pro Arg Val Gly Lys Ser Thr Ala Ile Leu Leu Gly Ile Thr Thr Leu Phe Thr Leu Leu Trp Ala Pro Arg Val Phe Val 280 Met Leu Tyr His Met Tyr Val Ala Pro Val His Arg Asp Trp Arg Val

His Leu Ala Leu Asp Val Ala Asn Met Val Ala Met Leu His Thr Ala

Ala Asn Phe Gly Leu Tyr Cys Phe Val Ser Lys Thr Phe Arg Ala Thr 325 330 335

Val Arg Gln Val Ile His Asp Ala Tyr Leu Pro Cys Thr Leu Ala Ser 340 345 350

Gln Pro Glu Gly Met Ala Ala Lys Pro Val Met Glu Pro Pro Gly Leu 355 360 365 <213>

Artificial Sequence

<220> <223> Novel Sequence <400> 7 atgtccactg aatgcgcgcg ggcagcgggc gacgcgccct tgcgcagcct ggagcaagcc 60 aaccgcaccc gctttccctt cttctccgac gtcaagggcg accaccggct ggtgctggcc 120 gcggtggaga caaccgtgct ggtgctcatc tttgcagtgt cgctgctggg caacgtgtgc 180 gccctggtgc tggtggcgcg ccgacgacgc cgcggcgcga ctgcctqcct gqtactcaac 240 ctettetgeg eggacetget etteateage getatecete tggtgetgge egtgegetgg 300 actgaggect ggctgctggg ccccgttgcc tgccacctgc tcttctacgt gatgaccctg 360 ageggcageg teaceateet eacgetggee geggteagee tggagegeat ggtgtgeate 420 gtgcacctgc agcgcggcgt gcggggtcct gggcggcggg cgcgggcagt gctgctggcg 480 ctcatatggg gctattcggc ggtcgccgct ctgcctctat gcgtcttctt ccgagtcgtc 540 ccgcaacggc tccccggcgc cgaccaggaa atttcgattt gcacactgat ttggcccacc 600 attcctggag agatctcgtg ggatgtctct tttgttactt tgaacttctt ggtgccagga 660 ctggtcattg tgatcagtta ctccaaaatt ttacagatca caaaggcatc aaggaagagg 720 ctcacggtaa gcctggccta ctcggagagc caccagatcc gcgtgtccca gcaggacttc 780 eggetettee geaccetett cetecteatg gteteettet teateatgtg gagececate 840 atcatcacca tectectcat cetgatecag aactteaage aagacetggt catetggeeg 900 tecetettet tetgggtggt ggeetteaea tttgetaatt cagecetaaa eeceateete 960 tacaacatga cactgtgcag gaatgagtgg aagaaaattt tttgctqctt ctgqttccca 1020

<210> 8 <211> 361 <212> PRT <213> Artificial Sequence <220>

<223> Novel Sequence

<400> 8

ggctaa

Met Ser Thr Glu Cys Ala Arg Ala Ala Gly Asp Ala Pro Leu Arg Ser 1 10 15

gaaaagggag ccattttaac agacacatct gtcaaaagaa atgacttgtc gattatttct

1080

1086

Leu Glu Gln Ala Asn Arg Thr Arg Phe Pro Phe Phe Ser Asp Val Lys Gly Asp His Arg Leu Val Leu Ala Ala Val Glu Thr Thr Val Leu Val Leu Ile Phe Ala Val Ser Leu Leu Gly Asn Val Cys Ala Leu Val Leu Val Ala Arg Arg Arg Arg Gly Ala Thr Ala Cys Leu Val Leu Asn 65 70 75 80 Leu Phe Cys Ala Asp Leu Leu Phe Ile Ser Ala Ile Pro Leu Val Leu Ala Val Arg Trp Thr Glu Ala Trp Leu Leu Gly Pro Val Ala Cys His Leu Leu Phe Tyr Val Met Thr Leu Ser Gly Ser Val Thr Ile Leu Thr Leu Ala Ala Val Ser Leu Glu Arg Met Val Cys Ile Val His Leu Gln Arg Gly Val Arg Gly Pro Gly Arg Arg Ala Arg Ala Val Leu Leu Ala Leu Ile Trp Gly Tyr Ser Ala Val Ala Ala Leu Pro Leu Cys Val Phe Phe Arg Val Val Pro Gln Arg Leu Pro Gly Ala Asp Gln Glu Ile Ser Ile Cys Thr Leu Ile Trp Pro Thr Ile Pro Gly Glu Ile Ser Trp Asp Val Ser Phe Val Thr Leu Asn Phe Leu Val Pro Gly Leu Val Ile Val Ile Ser Tyr Ser Lys Ile Leu Gln Ile Thr Lys Ala Ser Arg Lys Arg Leu Thr Val Ser Leu Ala Tyr Ser Glu Ser His Gln Ile Arg Val Ser Gln Gln Asp Phe Arg Leu Phe Arg Thr Leu Phe Leu Leu Met Val Ser 260 Phe Phe Ile Met Trp Ser Pro Ile Ile Ile Thr Ile Leu Leu Ile Leu Page 9

275 280 285

Ile Gln Asn Phe Lys Gln Asp Leu Val Ile Trp Pro Ser Leu Phe Phe 290 295 300

Trp Val Val Ala Phe Thr Phe Ala Asn Ser Ala Leu Asn Pro Ile Leu 305 310 315 320

Tyr Asn Met Thr Leu Cys Arg Asn Glu Trp Lys Lys Ile Phe Cys Cys 325 330 335

Phe Trp Phe Pro Glu Lys Gly Ala Ile Leu Thr Asp Thr Ser Val Lys 340 345 350

Arg Asn Asp Leu Ser Ile Ile Ser Gly

<210> 9

<211> 1038

<212> DNA

<213> Artificial Sequence

<220>

<223> Novel Sequence

<400> 9

atgagcagca attcatccct gctggtggct gtgcagctgt gctacgcgaa cgtgaatggg 60 tectgtgtga aaateeeett etegeeggga teeegggtga ttetgtacat agtgtttgge 120 tttggggctg tgctggctgt gtttggaaac ctcctggtga tgatttcaat cctccatttc 180 aagcagctgc actctccgac caattttctc gttgcctctc tggcctgcgc tgatttcttg 240 gtgggtgtga ctgtgatgcc cttcagcatg gtcaggacgg tggagagctg ctggtatttt 300 gggaggagtt tttgtacttt ccacacctgc tgtgatgtgg cattttgtta ctcttctctc 360 tttcacttgt gcttcatctc catcgacagg tacattgcgg ttactgaccc cctggtctat 420 cctaccaagt tcaccgtatc tgtgtcagga atttgcatca gcgtgtcctg gatcctgccc 480 ctcatgtaca gcggtgctgt gttctacaca ggtgtctatg acgatgggct ggaggaatta 540 tctgatgccc taaactgtat aggaggttgt cagaccgttg taaatcaaaa ctgggtgttg 600 acagattttc tatccttctt tatacctacc tttattatga taattctgta tggtaacata 660 tttcttgtgg ctagacgaca ggcgaaaaag atagaaaata ctggtagcaa gacagaatca 720 tcctcagaga gttacaaagc cagagtggcc aggagagag gaaaagcagc taaaaccctg 780 ggggtcacag tggtagcatt tatgatttca tggttaccat atagcattga ttcattaatt 840 gatgccttta tgggctttat aacccctgcc tgtatttatg agatttgctg ttggtgtqct 900 tattataact cagccatgaa tcctttgatt tatgctttat tttacccatg gtttaggaaa 960 qcaataaaag ttattgtaac tggtcaggtt ttaaagaaca gttcagcaac catgaatttg 1020

<210> 10

<211> 345

<212> PRT

<213> Artificial Sequence

<220>

<223> Novel Sequence

<400> 10

Met Ser Ser Asn Ser Ser Leu Leu Val Ala Val Gln Leu Cys Tyr Ala 1 5 10 15

Asn Val Asn Gly Ser Cys Val Lys Ile Pro Phe Ser Pro Gly Ser Arg 20 25 30

Val Ile Leu Tyr Ile Val Phe Gly Phe Gly Ala Val Leu Ala Val Phe 35 40 45

Gly Asn Leu Leu Val Met Ile Ser Ile Leu His Phe Lys Gln Leu His 50 55 60

Ser Pro Thr Asn Phe Leu Val Ala Ser Leu Ala Cys Ala Asp Phe Leu 65 70 75 80

Val Gly Val Thr Val Met Pro Phe Ser Met Val Arg Thr Val Glu Ser 85 90 95

Cys Trp Tyr Phe Gly Arg Ser Phe Cys Thr Phe His Thr Cys Cys Asp 100 105 110

Val Ala Phe Cys Tyr Ser Ser Leu Phe His Leu Cys Phe Ile Ser Ile 115 120 125

Asp Arg Tyr Ile Ala Val Thr Asp Pro Leu Val Tyr Pro Thr Lys Phe 130 135 140

Thr Val Ser Val Ser Gly Ile Cys Ile Ser Val Ser Trp Ile Leu Pro 145 150 155 160

Leu Met Tyr Ser Gly Ala Val Phe Tyr Thr Gly Val Tyr Asp Asp Gly
165 170 175

Leu Glu Glu Leu Ser Asp Ala Leu Asn Cys Ile Gly Gly Cys Gln Thr 180 185 190

Val Val Asn Gln Asn Trp Val Leu Thr Asp Phe Leu Ser Phe Phe Ile 195 200 205

Pro Thr Phe Ile Met Ile Ile Leu Tyr Gly Asn Ile Phe Leu Val Ala Page 11 Arg Arg Gln Ala Lys Lys Ile Glu Asn Thr Gly Ser Lys Thr Glu Ser 240

Ser Ser Glu Ser Tyr Lys Ala Arg Val Ala Arg Arg Glu Arg Lys Ala 250

Ala Lys Thr Leu Gly Val Thr Val Val Ala Phe Met Ile Ser Trp Leu 260 265 270

Pro Tyr Ser Ile Asp Ser Leu Ile Asp Ala Phe Met Gly Phe Ile Thr 275 280 285

Pro Ala Cys Ile Tyr Glu Ile Cys Cys Trp Cys Ala Tyr Tyr Asn Ser 290 295 300

Ala Met Asn Pro Leu Ile Tyr Ala Leu Phe Tyr Pro Trp Phe Arg Lys 305 310 315 320

Ala Ile Lys Val Ile Val Thr Gly Gln Val Leu Lys Asn Ser Ser Ala 325 330 335

Thr Met Asn Leu Phe Ser Glu His Ile 340 345

340 343

<210> 11

<211> 1020

<212> DNA

<213> Artificial Sequence

<220>

<223> Novel Sequence

<400> 11

atgatgccct tttgccacaa tataattaat atttcctgtg tgaaaaacaa ctggtcaaat 60 gatgtccgtg cttccctgta cagtttaatg gtgctcataa ttctgaccac actcgttggc 120 aatctgatag ttattgtttc tatatcacac ttcaaacaac ttcatacccc aacaaattgg 180 ctcattcatt ccatggccac tgtggacttt cttctggggt gtctggtcat gccttacagt 240 atggtgagat ctgctgagca ctgttggtat tttggagaag tcttctgtaa aattcacaca 300 agcaccgaca ttatgctgag ctcagcctcc attttccatt tgtctttcat ctccattgac 360 cgctactatg ctgtgtgtga tccactgaga tataaagcca agatgaatat cttggttatt 420 tgtgtgatga tcttcattag ttggagtgtc cctgctgttt ttgcatttgg aatgatcttt 480 ctggagctaa acttcaaagg cgctgaagag atatattaca aacatgttca ctgcagagga 540 ggttgctctg tcttctttag caaaatatct ggggtactga cctttatgac ttctttttat 600 atacctggat ctattatgtt atgtgtctat tacagaatat atcttatcgc taaagaacag 660

gca	agat	taa	ttag	tgat	gc c	aatc	agaa	g ct	ccaa	attg	gat	tgga	aat	gaaa	aatgga
att	tcac	aaa	gcaa	agaa	ag g	aaag	ctgt	g aa	gaca	ttgg	gga	ttgt	gat	ggga	gttttc
cta	atat	gct	ggtg	ccct	tt c	ttta	tctg	t ac	agtc	atgg	acc	cttt	tct	tcac	tacatt
att	ccac	cta	cttt	gaat	ga t	gtat	tgat	t tg	gttt	ggct	act	tgaa	ctc	taca	tttaat
cca	atgg	ttt i	atgc	attt	tt c	tatc	cttg	g tt	taga	aaag	cac	tgaa	gat	gatg	ctgttt
ggt	aaaa	ttt ·	tcca	aaaa	ga t	tcat	ccag	g tg	taaa	ttat	ttt	tgga	att	gagt	tcatag
<21 <21	<210> 12 <211> 339 <212> PRT <213> Artificial Sequence														
<22 <22		Nove:	l Se	quen	ce										
<40	0>	12							•						
Met 1	Met	Pro	Phe	Cys 5	His	Asn	Ile	Ile	Asn 10	Ile	Ser	Cys	Val	Lys 15	Asn
Asn	Trp	Ser	Asn 20	Asp	Val	Arg	Ala	Ser 25	Leu	Tyr	Ser	Leu	Met 30	Val	Leu
Ile	Ile	Leu 35	Thr	Thr	Leu	Val	Gly 40	Asn	Leu	Ile	Val	Ile 45	Va1	Ser	Ile
Ser	His 50	Phe	Lys	Gln	Leu	His 55	Thr	Pro	Thr	Asn	Trp 60	Leu	Ile	His	Ser
Met 65	Ala	Thr	Val	Asp	Phe 70	Leu	Leu	Gly	Суз	Leu 75	Val	Met	Pro	Tyr	Ser 80
Met	Val	Arg	Ser	Ala 85	Glu	His	Суз	Trp	Туг 90	Phe	Gly	Glu	Val	Phe 95	Cys
Lys	Ile	His	Thr 100	Ser	Thr	Asp	Ile	Met 105	Leu	Ser	Ser	Ala	Ser 110	Ile	Phe
His	Leu	Ser 115	Phe	Ile	Ser	Ile	Asp 120	Arg	Tyr	Tyr	Ala	Val 125	Cys	Asp	Pro
Leu	Arg 130	Tyr	Lys	Ala	Lys	Met 135	Asn	Ile	Leu	Val	Ile 140	Cys	Val	Met	Ile
Phe 145	Ile	Ser	Trp	Ser	Val 150	Pro	Ala	Val	Phe	Ala 155	Phe	Gly	Met	Ile	Phe 160

Leu Glu Leu Asn Phe Lys Gly Ala Glu Glu Ile Tyr Tyr Lys His Val 165 170 175

His Cys Arg Gly Gly Cys Ser Val Phe Phe Ser Lys Ile Ser Gly Val 180 185 190										
Leu Thr Phe Met Thr Ser Phe Tyr Ile Pro Gly Ser Ile Met Leu Cys 195 200 205										
Val Tyr Tyr Arg Ile Tyr Leu Ile Ala Lys Glu Gln Ala Arg Leu Ile 210 215 220										
Ser Asp Ala Asn Gln Lys Leu Gln Ile Gly Leu Glu Met Lys Asn Gly 235 240										
Ile Ser Gln Ser Lys Glu Arg Lys Ala Val Lys Thr Leu Gly Ile Val 245 250 255										
Met Gly Val Phe Leu Ile Cys Trp Cys Pro Phe Phe Ile Cys Thr Val 260 265 270										
Met Asp Pro Phe Leu His Tyr Ile Ile Pro Pro Thr Leu Asn Asp Val 275 280 285										
Leu Ile Trp Phe Gly Tyr Leu Asn Ser Thr Phe Asn Pro Met Val Tyr 290 295 300										
Ala Phe Phe Tyr Pro Trp Phe Arg Lys Ala Leu Lys Met Met Leu Phe 305 310 315 320										
Gly Lys Ile Phe Gln Lys Asp Ser Ser Arg Cys Lys Leu Phe Leu Glu 325 330 335										
Leu Ser Ser										
<210> 13 <211> 1029 <212> DNA <213> Artificial Sequence										
<220> <223> Novel Sequence										
<400> 13 atgaccagca atttttccca acctgttgtg cagctttgct atgaggatgt gaatggatct	60									
tgtattgaaa ctccctattc tcctgggtcc cgggtaattc tgtacacggc gtttagcttt 12	20									
gggtctttgc tggctgtatt tggaaatctc ttagtaatga cttctgttct tcattttaag 18	80									
	40									
	00									
gccaaatttt gtactettea eagttgetgt gatgtggeat tttgttacte ttetgteete 3	60									

cacttgtgct	tcatctgcat	cgacaggtac	attgtggtta	ctgatcccct	ggtctatgct	420
accaagttca	ccgtgtctgt	gtcgggaatt	tgcatcagcg	tgtcctggat	tctgcctctc	480
acgtacagcg	gtgctgtgtt	ctacacaggt	gtcaatgatg	atgggctgga	ggaattagta	540
agtgctctca	actgcgtagg	tggctgtcaa	attattgtaa	gtcaaggctg	ggtgttgata	600
gattttctgt	tattcttcat	acctaccctt	gttatgataa	ttctttacag	taagattttt	660
cttatagcta	aacaacaagc	tataaaaatt	gaaactacta	gtagcaaagt	agaatcatcc	720
tcagagagtt	ataaaatcag	agtggccaag	agagagagga	aagcagctaa	aaccctgggg	780
gtcacggtac	tagcatttgt	tatttcatgg	ttaccgtata	cagttgatat	attaattgat	840
gcctttatgg	gcttcctgac	ccctgcctat	atctatgaaa	tttgctgttg	gagtgcttat	900
tataactcag	ccatgaatcc	tttgatttat	gctctatttt	atccttggtt	taggaaagcc	960
ataaaactta	ttttaagtgg	agatgttta	aaggctagtt	catcaaccat	tagtttattt	1020
ttagaataa						1029
<210> 14						

<211> 342 <212> PRT

<213> Artificial Sequence

<220>

<223> Novel Sequence

<400> 14

Met Thr Ser Asn Phe Ser Gln Pro Val Val Gln Leu Cys Tyr Glu Asp

Val Asn Gly Ser Cys Ile Glu Thr Pro Tyr Ser Pro Gly Ser Arg Val

Ile Leu Tyr Thr Ala Phe Ser Phe Gly Ser Leu Leu Ala Val Phe Gly 35 40 45

Asn Leu Leu Val Met Thr Ser Val Leu His Phe Lys Gln Leu His Ser

Pro Thr Asn Phe Leu Ile.Ala Ser Leu Ala Cys Ala Asp Phe Leu Val 65 70 80

Gly Val Thr Val Met Leu Phe Ser Met Val Arg Thr Val Glu Ser Cys 85 90 95

Trp Tyr Phe Gly Ala Lys Phe Cys Thr Leu His Ser Cys Cys Asp Val

Ala Phe Cys Tyr Ser Ser Val Leu His Leu Cys Phe Ile Cys Ile Asp 115 120 125

Arg Tyr Ile Val Val Thr Asp Pro Leu Val Tyr Ala Thr Lys Phe Thr 135 Val Ser Val Ser Gly Ile Cys Ile Ser Val Ser Trp Ile Leu Pro Leu Thr Tyr Ser Gly Ala Val Phe Tyr Thr Gly Val Asn Asp Asp Gly Leu Glu Glu Leu Val Ser Ala Leu Asn Cys Val Gly Gly Cys Gln Ile Ile Val Ser Gln Gly Trp Val Leu Ile Asp Phe Leu Leu Phe Phe Ile Pro Thr Leu Val Met Ile Ile Leu Tyr Ser Lys Ile Phe Leu Ile Ala Lys Gln Gln Ala Ile Lys Ile Glu Thr Thr Ser Ser Lys Val Glu Ser Ser Ser Glu Ser Tyr Lys Ile Arg Val Ala Lys Arg Glu Arg Lys Ala Ala 245 250 255Lys Thr Leu Gly Val Thr Val Leu Ala Phe Val Ile Ser Trp Leu Pro Tyr Thr Val Asp Ile Leu Ile Asp Ala Phe Met Gly Phe Leu Thr Pro Ala Tyr Ile Tyr Glu Ile Cys Cys Trp Ser Ala Tyr Tyr Asn Ser Ala Met Asn Pro Leu Ile Tyr Ala Leu Phe Tyr Pro Trp Phe Arg Lys Ala Ile Lys Leu Ile Leu Ser Gly Asp Val Leu Lys Ala Ser Ser Ser Thr Ile Ser Leu Phe Leu Glu 340 <210> 15 <211> 1062 <212> DNA <213> Artificial Sequence <220> <223> Novel Sequence <400> 15

```
atggagcaca cgcacgccca cctcgcagcc aacagctcgc tqtcttqqtq qtcccccqqc
                                                                       60
tcggcctgcg gcttgggttt cgtgcccgtg gtctactaca gcctcttgct gtgcctcggt
                                                                      120
ttaccagcaa atatcttgac agtgatcatc ctctcccagc tggtggcaag aagacagaag
                                                                      180
tectectaca actatetett ggeacteget getgeegaca tettggteet ettitteata
                                                                      240
gtgtttgtgg acttcctgtt ggaagatttc atcttgaaca tgcagatgcc tcaggtcccc
                                                                      300
gacaagatca tagaagtgct ggaattctca tccatccaca cctccatatg gattactgta
                                                                      360
ccgttaacca ttgacaggta tatcgctgtc tgccacccgc tcaagtacca cacggtctca
                                                                      420
tacccagece geacceggaa agteattgta agtgtttaca teacctgett cetgaccage
                                                                      480
atcccctatt actggtggcc caacatctgg actgaagact acatcagcac ctctgtgcat
                                                                      540
cacgtcctca tctggatcca ctgcttcacc gtctacctgg tgccctgctc catcttcttc
                                                                      600
atcttgaact caatcattgt gtacaagctc aggaggaaga gcaattttcg tctccgtggc
                                                                      660
tactccacgg ggaagaccac cgccatcttg ttcaccatta cctccatctt tgccacactt
                                                                      720
tgggccccc gcatcatcat gattctttac cacctctatg gggcgcccat ccagaaccgc
                                                                      780
tggctggtgc acatcatgtc cgacattgcc aacatgctaq cccttctgaa cacagccatc
                                                                      840
aacttettee tetaetgett cateageaag eggtteegea ceatggeage egeeaegete
                                                                      900
aaggetttet teaagtgeea gaageaacet gtacagttet acaccaatea taacttttee
                                                                      960
ataacaagta gcccctggat ctcgccggca aactcacact gcatcaagat gctggtgtac
                                                                    1020
cagtatgaca aaaatggaaa acctataaaa gtatccccgt ga
                                                                    1062
```

<210> 16

<211> 353

<212> PRT

<213> Artificial Sequence

<220>

<223> Novel Sequence

<400> 16

Met Glu His Thr His Ala His Leu Ala Ala Asn Ser Ser Leu Ser Trp 1 5 10 15

Trp Ser Pro Gly Ser Ala Cys Gly Leu Gly Phe Val Pro Val Val Tyr 20 25 30

Tyr, Ser Leu Leu Cys Leu Gly Leu Pro Ala Asn Ile Leu Thr Val 35 40 45 .

Ile Ile Leu Ser Gln Leu Val Ala Arg Arg Gln Lys Ser Ser Tyr Asn 50 60

Tyr Leu Leu Ala Leu Ala Ala Ala Asp Ile Leu Val Leu Phe Phe Ile 65 70 75 80

Val Phe Val Asp Phe Leu Leu Glu Asp Phe Ile Leu Asn Met Gln Met Pro Gln Val Pro Asp Lys Ile Ile Glu Val Leu Glu Phe Ser Ser Ile His Thr Ser Ile Trp Ile Thr Val Pro Leu Thr Ile Asp Arg Tyr Ile Ala Val Cys His Pro Leu Lys Tyr His Thr Val Ser Tyr Pro Ala Arg Thr Arg Lys Val Ile Val Ser Val Tyr Ile Thr Cys Phe Leu Thr Ser Ile Pro Tyr Tyr Trp Pro Asn Ile Trp Thr Glu Asp Tyr Ile Ser Thr Ser Val His His Val Leu Ile Trp Ile His Cys Phe Thr Val Tyr Leu Val Pro Cys Ser Ile Phe Phe Ile Leu Asn Ser Ile Ile Val Tyr Lys Leu Arg Arg Lys Ser Asn Phe Arg Leu Arg Gly Tyr Ser Thr Gly Lys Thr Thr Ala Ile Leu Phe Thr Ile Thr Ser Ile Phe Ala Thr Leu Trp Ala Pro Arg Ile Ile Met Ile Leu Tyr His Leu Tyr Gly Ala Pro Ile Gln Asn Arg Trp Leu Val His Ile Met Ser Asp Ile Ala Asn Met Leu Ala Leu Leu Asn Thr Ala Ile Asn Phe Phe Leu Tyr Cys Phe Ile 280 Ser Lys Arg Phe Arg Thr Met Ala Ala Ala Thr Leu Lys Ala Phe Phe Lys Cys Gln Lys Gln Pro Val Gln Phe Tyr Thr Asn His Asn Phe Ser Ile Thr Ser Ser Pro Trp Ile Ser Pro Ala Asn Ser His Cys Ile Lys 330 Met Leu Val Tyr Gln Tyr Asp Lys Asn Gly Lys Pro Ile Lys Val Ser Page 18

340 345 350

Pro

<210> 17 <211> 969 <212> DNA <213> Artificial Sequence									
<220> <223> Novel Sequence									
<400> 17 atggatccaa ccgtcccagt cttcggtaca aaactgacac caatcaacgg acgtgaggag	60								
actccttgct acaatcagac cctgagcttc acggtgctga cgtgcatcat ttcccttgtc	120								
ggactgacag gaaacgcggt agtgctctgg ctcctgggct accgcatgcg caggaacgct	180								
gtctccatct acatcctcaa cctggccgca gcagacttcc tcttcctcag cttccagatt	240								
atacgttcgc cattacgcct catcaatatc agccatctca tccgcaaaat cctcgtttct	300								
gtgatgacct ttccctactt tacaggcctg agtatgctga gcgccatcag caccgagcgc	360								
tgcctgtctg ttctgtggcc catctggtac cgctgccgcc gccccacaca cctgtcagcg	420								
gtcgtgtgtg tcctgctctg gggcctgtcc ctgctgttta gtatgctgga gtggaggttc	480								
tgtgacttcc tgtttagtgg tgctgattct agttggtgtg aaacgtcaga tttcatccca	540								
gtcgcgtggc tgatttttt atgtgtggtt ctctgtgttt ccagcctggt cctgctggtc	600								
aggatectet gtggateceg gaagatgeeg etgaceagge tgtacgtgae cateetgete	660								
acagtgctgg tcttcctcct ctgcggcctg cccttcggca ttctgggggc cctaatttac	720								
aggatgcacc tgaatttgga agtcttatat tgtcatgttt atctggtttg catgtccctg	780								
tcctctctaa acagtagtgc caaccccatc atttacttct tcgtgggctc ctttaggcag	840								
cgtcaaaata ggcagaacct gaagctggtt ctccagaggg ctctgcagga caagcctgag	900								
gtggataaag gtgaagggca gcttcctgag gaaagcctgg agctgtcggg aagcagattg	960								
gggccatga	969								
<210> 18 <211> 322 <212> PRT <213> Artificial Sequence									
<220> <223> Novel Sequence									
<400> 18									

Gly Arg Glu Glu Thr Pro Cys Tyr Asn Gln Thr Leu Ser Phe Thr Val Page 19 $\,$

Met Asp Pro Thr Val Pro Val Phe Gly Thr Lys Leu Thr Pro Ile Asn 1 5

- Leu Thr Cys Ile Ile Ser Leu Val Gly Leu Thr Gly Asn Ala Val Val 35 40 45
- Leu Trp Leu Leu Gly Tyr Arg Met Arg Arg Asn Ala Val Ser Ile Tyr 50 60
- Ile Leu Asn Leu Ala Ala Ala Asp Phe Leu Phe Leu Ser Phe Gln Ile 65 70 75 80
- Ile Arg Ser Pro Leu Arg Leu Ile Asn Ile Ser His Leu Ile Arg Lys 85 90 95
- Ile Leu Val Ser Val Met Thr Phe Pro Tyr Phe Thr Gly Leu Ser Met 100 105 110
- Leu Ser Ala Ile Ser Thr Glu Arg Cys Leu Ser Val Leu Trp Pro Ile 115 120 · 125
- Trp Tyr Arg Cys Arg Arg Pro Thr His Leu Ser Ala Val Val Cys Val 130 135 140
- Leu Leu Trp Gly Leu Ser Leu Leu Phe Ser Met Leu Glu Trp Arg Phe 145 150 160
- Cys Asp Phe Leu Phe Ser Gly Ala Asp Ser Ser Trp Cys Glu Thr Ser 165 170 175
- Asp Phe Ile Pro Val Ala Trp Leu Ile Phe Leu Cys Val Val Leu Cys 180 185 190
- Val Ser Ser Leu Val Leu Leu Val Arg Ile Leu Cys Gly Ser Arg Lys
- Met Pro Leu Thr Arg Leu Tyr Val Thr Ile Leu Leu Thr Val Leu Val 210 215 220
- Phe Leu Cys Gly Leu Pro Phe Gly Ile Leu Gly Ala Leu Ile Tyr 225 230 235 240
- Arg Met His Leu Asn Leu Glu Val Leu Tyr Cys His Val Tyr Leu Val
 245 250 255
- Cys Met Ser Leu Ser Ser Leu Asn Ser Ser Ala Asn Pro Ile Ile Tyr 260 265 270
- Phe Phe Val Gly Ser Phe Arg Gln Arg Gln Asn Arg Gln Asn Leu Lys 275 280 285

Leu Val Leu Gln Arg Ala Leu Gln Asp Lys Pro Glu Val Asp Lys Gly 290 Glu Gly Gln Leu Pro Glu Glu Ser Leu Glu Leu Ser Gly Ser Arg Leu 305 310 315 Gly Pro <210> 19 <211> 969 <212> <213> Artificial Sequence <220> <223> Novel Sequence <400> 19 atggattcaa ccatcccagt cttgggtaca gaactgacac caatcaacgg acgtgaggag 60 actecttget acaagcagae eetgagette acggggetga egtgcategt tteeettgtb 120 gcgctgacag gaaacgcggt tgtgctctgg ctcctgggct gccgcatgcg caggaacgct 180 gtctccatct acatcctcaa cctggtcgcg gccgacttcc tcttccttag cggccacatt 240 atatgttcgc cgttacgcct catcaatatc cgccatccca tctccaaaat cctcagtcct 300 gtgatgacct ttccctactt tataggccta agcatgctga gcgccatcag caccgagcgc 360 tgcctgtcca tcctgtggcc catctggtac cactgccgcc gccccagata cctgtcatcg 420 gtcatgtgtg tcctgctctg ggccctgtcc ctgctgcgga gtatcctgga gtggatgttc 480 tgtgacttcc tgtttagtgg tgctgattct gtttggtgtg aaacgtcaga tttcattaca 540 atcgcgtggc tggttttttt atgtgtggtt ctctgtgggt ccagcctggt cctgctggtc 600 aggattetet gtggateceg gaagatgeeg etgaceagge tgtacgtgae eatecteete 660 acagtgctgg tettectect etgtggcetg ceetttggca tteagtggge eetgttttee 720 aggatccacc tggattggaa agtcttattt tgtcatgtgc atctagtttc cattttcctg 780 tecgetetta acageagtge caaceceate atttacttet tegtgggete etttaggeag 840 cgtcaaaata ggcagaacct gaagctggtt ctccagaggg ctctgcagga cacgcctgag 900 gtggatgaag gtggagggtg gcttcctcag gaaaccctgg agctgtcggg aagcagattg 960 gagcagtga 969 <210> 20 <211> 322 <212> PRT <213> Artificial Sequence

<220>

<400> 20

<223> Novel Sequence

Met 1	Asp	Ser	Thr	Ile 5	Pro	Val	Leu	Gly	Thr 10	Glu	Leu	Thr	Pro	Ile 15	Asn
Gly	Arg	Glu	Glu 20	Thr	Pro	Суз	Tyr	Lys 25	Gln	Thr	Leu	Ser	Phe 30	Thr	Gly
Leu	Thr	Суз 35	Ile	Val	Ser	Leu	Val 40	Ala	Leu	Thr	Gly	Asn 45	Ala	Val	Val
Leu	Trp 50	Leu	Leu	Gly	Суз	Arg 55	Met	Arg	Arg	Asn	Ala 60	Val	Ser	Ile	Tyr
Ile 65	Leu	Asn	Leu	Val	Ala 70	Ala	Asp	Phe	Leu	Phe 75	Leu	Ser	Gly	His	Ile 80
Ile	Cys	Ser	Pro	Leu 85	Arg	Leu	Ile	Asn	Ile 90	Arg	His	Pro	Ile	Ser 95	Lys
Ile	Leu	Ser	Pro 100	Val	Met	Thr	Phe	Pro 105	Tyr	Phe	Ile	Gly	Leu 110	Ser	Met
Leu	Ser	Ala 115	Ile	Ser	Thr	Glu	Arg 120	Суз	Leu	Ser	Ile	Leu 125	Trp	Pro	Ile
Trp	Tyr 130	His	Cys	Arg	Arg	Pro 135	Arg	Tyr	Leu	Ser	Ser 140	Val	Met	Cys	Val
Leu 145	Leu	Trp	Ala	Leu	Ser 150	Leu	Leu	Arg	Ser	Ile 155	Leu	Glu	Trp	Met	Phe 160
Суѕ	Asp	Phe	Leu	Phe 165	Ser	Gly	Ala	Asp	Ser 170	Val	Trp	Cys	Glu	Thr 175	Ser
Asp	Phe	Ile	Thr 180	Ile	Ala	Trp	Leu	Val 185	Phe	Leu	Cys	Val	Val 190	Leu	Cys
Gly	Ser	<i>S</i> er 195	Leu	Val	Leu	Leu	Val 200	Arg	Ile	Leu	Cys	Gly 205	Ser	Arg	Lys
Met	Pro 210	Leu	Thr	Arg	Leu	Tyr 215	Val	Thr	Ile	Leu	Leu 220	Thr	Val	Leu	Val
Phe 225	Leu	Leu	Cys	Gly	Leu 230	Pro	Phe	Gly	Ile	Gln 235	Trp	Ala	Leu	Phe	Ser 240
Arg	Ile	His	Leu	Asp 245	Trp	Lys	Val	Leu	Phe 250	Cys	His	Val	His	Leu 255	Val
Ser	Ile	Phe	Leu	Ser	Ala	Leu	Asn	Ser	Ser		Asn Page		Ile	Ile	Tyr

Phe Ph	e Val 6 275	Gly Ser	Phe	Arg	Gln 280	Arg	Gln	Asn	Arg	Gln 285	Asn	Leu	Lys	
Leu Va 29		Gln Arg	Ala	Leu 295	Gln	Asp	Thr	Pro	Glu 300	Val	Asp	Glu	Gly	
Gly Gl 305	y Trp L	Leu Pro	Gln 310	Glu	Thr	Leu	Glu	Leu 315	Ser	Gly	Ser	Arg	Leu 320	
Glu Gl	n													
<210> <211> <212> <213>	21 26 DNA Artifi	icial S	equer	ıce										
<220> <223>	Novel	Sequen	ce											
<400> cagage	21 tctg gt	ggccac	ct ct	gtco	:									26·
<210> <211> <212> <213>		icial Se	equer	ıce									•	
<220> <223>	Novel	Sequenc	ce											
<400> ctgcgt	22 ccac ca	agagtca	cg to	tcc										25
<210> <211> <212> <213>	23 26 DNA Artifi	icial Se	equen	ıce										
<220> <223>	Novel	Sequenc	ce											
<400> gtatgc	23 ctgg cc	cacaatao	ec to	cagg	•									26
<210> <211> <212> <213>	24 31 DNA Artifi	icial Se	equen	ıce										
<220> <223>	Novel	Sequenc	ce ·											
<400> gtttgt	24 ggct aa	cggcaca	aa aa	caca	atto	c c								31

<210> <211> <212> <213>		
<220>		
<223>	Novel Sequence	
<400> ggtacc	25 acaa tgacaatcac cagcgtcc	28
<210> <211>	26 29	
<212> <213>	DNA	
	Artificial Sequence	
<220> <223>	Novel Sequence	
<400>		
ggaacg	tgag gtacatgtgg atgtgcagc	29
<210>	27	
<211> <212>		
<213>	Artificial Sequence	
<220> <223>	Novel Sequence	
<400> gcagtg	27 tagc ggtcaaccgt gagcagg	27
<210>	28	
<211> <212>	32 DNA	
<213>	Artificial Sequence	
<220> <223>	Novel Sequence	
<400>	28	20
cyayca	ggat ggcgatccag actgaggcgt gg	32
<210>	29	
<211> <212>	25 DNA	
<213>	Artificial Sequence	
<220> <223>	Novel Sequence	
<400>	29 cagc tggcgatgct gacag	25
3~39000		23
<210>	30	
<211> <212>	25 DNA	
<213>	Artificial Sequence	

<220> <223>	Novel Sequence	
<400>	30	0.5
grygee	atga gccaccctga gctcc	25
<210> <211>	31 24	
<212>		
<213>	Artificial Sequence	
<220>		
<223>	Novel Sequence	
<400>	31	0.4
gyaary	tcca ctgaatgege gegg	24
<210>	32	
<211> <212>	25 DNA	
<213>	Artificial Sequence	
<220>		
<223>	Novel Sequence	
<400>	32	
agctcg	ccag gtgtgagaaa ctcgg	25
<210>	33	
<211>	30	
<212> <213>	DNA Artificial Sequence	
12107	naturatur boquenos	
<220>	Name I Commence	
<223>	Novel Sequence	
<400>	33	
gcgtta	cgag cagcaattca tecetgetgg	30
<210>	34	
<211>	28	
	DNA	
<213>	Artificial Sequence	
<220> <223>	Novel Sequence	
	34	
<400> gtatcct	gaa cttcgtctat acaactgc	28
<210>	35	
<211>	29 DNA	
<212> <213>	DNA Artificial Sequence	
<220>		
<223>	Novel Sequence	
<400>	35	
ccctcac	gaa tgatgccctt ttgccacaa	29

```
<210> 36
<211> 28
<212> DNA
<213> Artificial Sequence
<220>
<223> Novel Sequence
<400> 36
atccatgtgg ttggtgcatg tggttcgt
                                                                     28
<210> 37
<211> 29
<212> DNA
<213> Artificial Sequence
<220>
<223> Novel Sequence
<400> 37
aaacaacaaa cagcagaacc atgaccagc
                                                                     29
<210> 38
<211> 30
<212> DNA
<213> Artificial Sequence
<220>
<223> Novel Sequence
<400> 38
                                                                     30
acatagagac aagtgacatg tgtgaaccac
<210> 39
<211> 26
<212> DNA
<213> Artificial Sequence
<220>
<223> Novel Sequence
<400> 39
ggtatgagac cgtgtggtac ttgagc ·
                                                                     26
<210> 40
<211> 29
<212> DNA
<213> Artificial Sequence
<220>
<223> Novel Sequence
<400> 40
gtggcagaca gcgatatacc tgtcaatgg
                                                                     29
<210> 41
<211>
      27
<212> DNA
<213> Artificial Sequence
```

<220> <223>	Novel Sequence	
<400> gcgcto	41 catgg agcacacgca cgcccac	27
<210> <211> <212> <213>	25 DNA	
<220> <223>	Novel Sequence	
<400> gaggca	42 gtag ttgccacacc tatgg	25
<210> <211> <212> <213>	29 DNA	
<220> <223>	Novel Sequence	
<400> catctg	43 gttt gtgttcccag gggcaccag	29
<210> <211> <212> <213>	32	
<220> <223>	Novel Sequence	
	44 gttg ctctcaaagt cccgtctgac tg	32
<210> <211> <212> <213>	45 25 DNA Artificial Sequence	
<220> <223>	Novel Sequence	
<400> ctgttt	45 ccag ggtcatcaga ctggg	25
<210> <211> <212> <213>	46 27 DNA Artificial Sequence	
<220> <223>	Novel Sequence	
<400>	46	27

```
<210> 47
<211> 6
<212> PRT
<213> Artificial Sequence
<220>
<223> Novel Sequence
<400> 47
Thr Leu Glu Ser Ile Met
<210> 48
<211> 5
<212> PRT
<213> Artificial Sequence
<220>
<223> Novel Sequence
<400> 48
Glu Tyr Asn Leu Val
<210> 49
<211> 5
<212> PRT
<213> Artificial Sequence
<220>
<223> Novel Sequence
<400> 49
Asp Cys Gly Leu Phe
<210> 50
<211> 34
<212> DNA
<213> Artificial Sequence
<220>
<223> Novel Sequence
<400> 50
                                                                         34
gatcaagctt ccatggcgtg ctgcctgagc gagg
<210> 51
<211> 53
<212> DNA
<213> Artificial Sequence
<220>
<223> Novel Sequence
<400> 51
                                                                         53
gateggatec ttagaacagg cegeagtect teaggtteag etgeaggatg gtg
```

<210> <211> <212> <213>	52 25 DNA Artificial Sequence	
<220> <223>	Novel Sequence	
<400> gtcctc	52 actg gtggccatgt actcc	25
<210> <211> <212> <213>		
<220> <223>	Novel Sequence	
	53 ccac cagagtcacg tetec	25
<210> <211> <212> <213>		
<220> <223>	Novel Sequence	
	54 tgtt tgggctgccc ttctgc	26
<210> <211> <212> <213>	DNA	
<220> <223>	Novel Sequence	
<400> gtttgt	55 gget aacggcacaa aacacaatte c	31
<210> <211> <212> <213>	56 26 DNA Artificial Sequence	
<220> <223>	Novel Sequence	
<400> ctgctc	56 acgg ttgaccgcta cactgc	26
<210> <211> <212> <213>		

	<220> <223>	Novel Sequence	
	<400> stggcc	57 atga gccaccctga gctcc	25
<	<210> 1 1 1 1 1 1 1 1 1 1</td <td>20</td> <td></td>	20	
	(220> (223>	Novel Sequence	
	(400> ettett	58 . ctcc gacgtcaagg	20
<	(210> (211> (212> (213> (220> (220>	20 DNA	
<	(400>		20
<	220>	20	
	400> ttctt	60 etcc gacgtcaagg	20
<	210> 211> 212> 213>	20 .	
	220> 223>	Novel Sequence	
	400> aacgg	61 totg acaacctcct	20
<	210> 211> 212> 213>	DNA	•
	220>	Novel Sequence	
	400> tgctg	62 tgat gtggcatttt g	21

<210> <211> <212> <213>	DNA	
<220> <223>	Novel Sequence	
<400> caggaa	63 gccc ataaaggcat caa	23
<210> <211> <212> <213>		
<220> <223>	Novel Sequence	
<400> acatca	64 cctg cttcctgacc	20
<210> <211> <212> <213>	65 20 DNA Artificial Sequence	
<220> <223>	Novel Sequence	
<400> ccagca	65 tctt gatgcagtgt	20
<210> <211> <212> <213>	66 21 DNA Artificial Sequence	
<220> <223>	Novel Sequence	
<400> ccatct	66 ccaa aatcctcagt c	21
<210> <211> <212> <213>	67 22 DNA Artificial Sequence	
<220> <223>	Novel Sequence	
<400> gctgtt	67 aaga gcggacagga aa	22